

# Making a Difference in the Journey of Life



## Ostomy Care Product Catalog



Ostomy Care  
*Healthy Skin. Positive Outcomes.*

# Welcome to your Hollister Ostomy Product Guide

Within this guide you will find a range of products designed to suit individual needs, including CeraPlus™ Ostomy Products\*, an advanced line of ostomy products that help keep skin healthy.



CeraPlus™ Ostomy Products have received official dermatological accreditation from the Skin Health Alliance, bringing confidence and reassurance to end users.

Healthy Skin. Positive Outcomes.

# Promoting Healthy Peristomal Skin

## Promoting healthy peristomal skin

It drives everything we do. In order to promote healthy skin, it is important to understand that the fit of the pouching system as well as the formulation of the skin barrier go hand in hand.

A pouching system will not maintain a good seal without the proper fit. Without the correct formulation, the skin around the stoma can be at risk for complications.

Our goal is to help keep the skin around your patient's stoma as healthy as the skin on the rest of their body.

<b>CeraPlus™ Products</b> .....	13-28
Two-Piece Skin Barriers .....	15-16
Two-Piece Drainable Pouches .....	18-20
Two-Piece Closed Pouches .....	19-20
Two-Piece Urostomy Pouches.....	20
One-Piece Drainable Pouches .....	23-26
One-Piece Urostomy Pouches .....	27-28
<b>Two-Piece Pouching Systems</b> .....	33-36
<b>FlexTend™</b> (Extended wear) Skin Barriers.....	34-35
<b>FlexWear™</b> (Standard wear) Skin Barriers.....	36
<b>One-Piece Pouching Systems</b> .....	37-42
<b>FlexTend™</b> (Extended wear) Skin Barriers.....	37-40
<b>FlexWear™</b> (Standard wear) Skin Barriers.....	41
<b>SoftFlex™</b> (Standard wear) Skin Barriers .....	42
<b>Hospital Items</b> .....	45-50
High Output.....	45-46
Wound Drainage.....	47
Sterile/Non-Sterile Kits .....	48-50
<b>Pediatrics</b> .....	51-54
<b>Accessories</b> .....	55-65
Ceramide Infused .....	56
Barrier Rings .....	56-58
Other Accessories.....	58-63
Irrigation Products .....	64-65
<b>Education</b> .....	69-73
<b>Index</b> .....	74-75
<b>HCPC / Medicare Utilization Chart</b> .....	76-81

# Hollister Secure Start™ Services

Personalized Support

## Why Secure Start™ Services for Your Patients?



**Secure Start™ services** helps patients—regardless of product brand they use—to get back to living life with confidence and security. We understand this is a process that takes time. That's why our no-cost, personalized service is available for the lifetime of their ostomy.



**THE RIGHT PRODUCTS**—based on your product recommendations, we help patients find the ostomy products that best fit their needs throughout the lifetime of their ostomy.



**OSTOMY SMART EDUCATION**—easy access to important care information and lifestyle resources that support the goal of maintaining healthy peristomal skin from the start.



**SUPPLIER MATCHING**—connects patients with a supplier of their choice based on their insurance plan for the products you recommend.



**WOC ACCESS**—if the need arises, an ostomy nurse is available to assist with product questions.

### Support Throughout the Continuum of Care

We are proud to offer dedicated support for each and every patient we serve along the continuum of care. Hollister Secure Start services provide a lifetime of personalized support for our members.

Our team is focused on helping patients with:

- Finding the right product fit
- Navigating insurance coverage
- Identifying supplier options
- Providing condition specific education
- And more!

To learn more, call us at **1.888.808.7456**.

### ENROLLMENT MADE EASY

Obtaining written patient consent to enroll in Secure Start™ services is an important step. At Hollister, we prioritize patient privacy—giving patients comfort that their information will remain confidential and will be shared only with those they wish to hear from. Having patient consent on file helps us provide proactive ostomy education and product information throughout their journey.

#### Digital

Easy mobile enrollment with Enroll on the Go



Scan Me

#### Fax

Fax enrollment form to 1.800.398.2341.

If you have any questions about the enrollment process, please call your Hollister Secure Start™ services coordinator at 1.888.808.7456. Learn more at [www.securestartservices.com](http://www.securestartservices.com).



Hollister, the Hollister logo and Secure Start are trademarks of Hollister Incorporated  
© 2024 Hollister Incorporated.




# CeraPlus™ Ostomy Products

Protection where it matters most



CeraPlus™ and barriers contain the Remois technology of Remois Co., Ltd.  
Please to use, be sure to read the instructions for use for information regarding Intended Use, Contraindications, Warnings, Precautions, and Instructions for Use. Hollister, the Hollister logo and CeraPlus are trademarks of Hollister Incorporated. All other trademarks and copyrights are the property of their respective owners. © 2024 Hollister Incorporated.

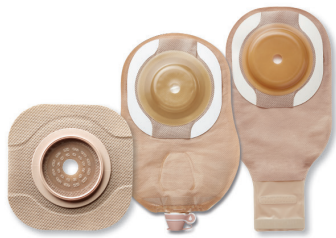

## Start with CeraPlus™ Skin Barriers\*

Protecting peristomal skin begins with a skin barrier that provides a secure, comfortable fit and a skin-friendly formulation. We focus on security and skin health when designing our skin barriers – including the ceramide-infused CeraPlus™ Portfolio.

The CeraPlus™ Line of Ostomy Products provide skin with advanced protection, and a full range of ceramide-infused options so you can create the personalized solution that's best for your patient's body.



Ostomy Care  
Healthy skin. Positive outcomes.



CeraPlus™ and barriers contain the Remois technology of Remois Co., Ltd.

**FIT** is critical to get a good seal around the stoma to prevent leakage


## Barrier Examples



**Flat Skin Barrier**



**Integrated Soft Convex Barrier**

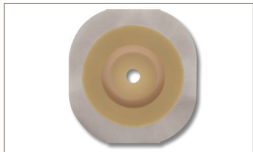

**Integrated Convex Barrier**



**Integrated Convex Barrier with Round or Oval Convex Barrier Ring**

*Round barrier ring shown in picture.*

## Possible Clinical Case Scenarios



**Flat Skin Barrier**

- Flat peristomal plane
- Adequate stoma protrusion
- Centered stoma opening




**Soft Convexity**

- Creases or folds in peristomal plane
- Gentle pressure required to peristomal plane
- Flush stoma
- Peristomal injury



**Firm Convexity**

- Creases or folds in peristomal plane
- Off-centered stoma opening
- Stoma retracted below skin surface

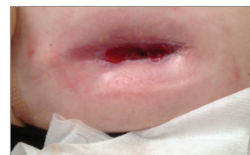

**Firm Convexity with More Depth**

- Recessed or flushed stoma
- Uneven abdominal contours
- Off-centered stoma opening

Prior to use, be sure to read the Instructions for Use for information regarding intended use, contraindications, warnings, precautions and instructions.

**FORMULATION** is the customized ingredients responsible for adhesion, erosion resistance and absorption that must be carefully formulated to create an environment where healthy skin can thrive.

## Types of Skin Barrier Formulations and Their Performance

The Hollister Ostomy portfolio has a multitude of formulation options that are designed to help ensure an optimal skin interface is possible for every patient.

### Extended Wear

#### CeraPlus™ Skin Barrier

- Contains all skin barrier characteristics (adhesion, absorption and erosion resistance) and is infused with ceramide – a naturally occurring lipid in the skin which protects against dryness
- Its formulation is designed to help protect the skin's natural moisture barrier and help maintain good peristomal skin health

#### FlexWear™ Skin Barrier

- Durable with high resistance to erosion
- Swells when it comes in contact with liquid to help form a seal between the skin and the stoma
- High initial tack and suitable for those who perspire heavily
- Strong adhesion and flexibility

#### SoftFlex™ Skin Barrier

### Standard Wear

#### FlexWear™ Skin Barrier

- Good balance of durability and gentleness providing moderate adhesion, and can handle frequent or infrequent pouch changes
- Resistant to fluid erosion and tends to hold its shape in the presence of fluids
- Suitable for all skin types, depending on change frequency
- Designed to be worn more than one day

#### SoftFlex™ Skin Barrier

- Appropriate choice for:
  - Those who have thin, fragile, or sensitive skin (such as those who have undergone chemo/radiation therapy)
  - People who desire frequent removal
- Designed to limit skin damage as it's gentle to the skin, which allows for frequent pouch changes
- Thin and flexible without compromising security

## CeraPlus™ Skin Barriers are demonstrated by evidence to reduce the occurrence and impact of PSCs<sup>1</sup>

As many as 75% of ostomates have peristomal skin complications (PSCs).<sup>1</sup> CeraPlus™ skin barriers protect skin from Day 1 making it an ideal first-line choice for your ostomy patients. When compared to non CeraPlus™ skin barriers, significantly more (75%) patients were “very satisfied” with the overall performance of CeraPlus™ skin barriers.<sup>1</sup>

- Helps support and maintain healthy peristomal skin<sup>1</sup>
- Delivers positive outcomes for ostomy-related costs and care<sup>1</sup>
- Provide fit with convexity options and adhesion performance to enhance skin health



Contact your sales specialist for free samples



1. Colwell J, Pittman J, Raizman R, Salvadalena G. A Randomized Controlled Trial Determining Variances in Ostomy Skin Conditions (ADVOCATE). J Wound Ostomy Continence Nurs. 2018;45(1):37-42.



## Security + Skin Health = Protection Where it Matters Most

Protecting peristomal skin means using products that provide a good seal around the stoma and promote skin health. Hollister uses the unique combination of security and skin health to design skin barriers that help maintain healthy skin – including the ceramide-infused CeraPlus™ Portfolio.\*





# CeraPlus™ Products

Infused with Ceramide



## Transforming the skin health conversation


At Hollister Ostomy Care, we believe peristomal skin deserves better. In fact, we think we should expect the skin around the stoma to be as healthy as the skin on the other side of the abdomen.

That's why we've made the pursuit of healthy skin our priority.

The CeraPlus™ Skin Barriers are formulated to support healthy skin from Day One. The ceramide-infused barriers are designed to help protect the skin's natural moisture barrier. Healthy skin can remain healthy through pouching changes, and support an appropriate fit around the stoma.



## **CeraPlus™ Skin Barriers - Available in One-Piece and Two-Piece Pouching Systems**

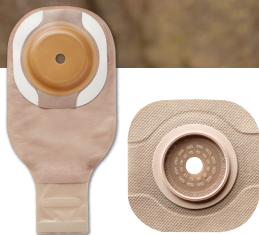

At Hollister, we realize the importance of healthy peristomal skin. Our CeraPlus™ skin barrier is infused with ceramide, the skin's naturally occurring protection against dryness. CeraPlus™ skin barrier is designed to help protect the skin's own moisture barrier, and is available in one-piece and two-piece pouching systems. This means you can enjoy the same pouch and skin barrier features from Hollister that you've come to rely on—and more!

*The CeraPlus™ Portfolio is the only advanced line of ostomy products in the U.S. infused with ceramide that protects skin from Day One.\**



\*CeraPlus Products contain the Remois technology of Alcare Co., LTD  
Hollister Data on File. Ref-02739; 2021.



**Ostomy Care**  
Healthy Skin. Positive Outcomes.

**CeraPlus™ Skin Barrier delivers both Security and Skin Health**

**CeraPlus™ Skin Barrier**

- Protects the skin's natural moisture barrier
- Decreases TEWL (transepidermal water loss) from damaged or eroded skin
- Maintains peristomal skin health





**Color Matching System**

**Two-Piece Pouching Systems**

To help you select compatible two-piece pouching system components, every box label is color-coded.

1. Select the skin barrier flange size with the appropriate barrier stoma opening and note its color code.
2. Select a suitable pouch with the desired features in the matching color.

**Flange Sizes:**

	1-3/4" (44 mm)
	2-1/4" (57 mm)
	2-3/4" (70 mm)
	4" (102 mm)

**Look for the convex and soft convex symbols:**





**CeraPlus™ Soft Convex Skin Barriers with Integrated Floating Flange**



**1170x (CUT-TO-FIT) Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Soft Convex
- ▶ Adhesive Border (with tape)

**HCPCS: A4407**





Stock No	Flange Size	For Stomas Up To
11702	 1-3/4" (44 mm)	1" (25 mm)
11703	 2-1/4" (57 mm)	1-1/2" (38 mm)
11704	 2-3/4" (70 mm)	2" (51 mm)



**1190x (PRE-SIZED) Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Soft Convex
- ▶ Adhesive Border (with tape)

**HCPCS: A4407**

Stock No	Flange Size	Barrier Opening
11903	 1-3/4" (44 mm)	7/8" (22 mm)
11904	 1-3/4" (44 mm)	1" (25 mm)
11905	 2-1/4" (57 mm)	1-1/8" (29 mm)
11906	 2-1/4" (57 mm)	1-1/4" (32 mm)




**CeraPlus™ Convex Skin Barriers with Integrated Floating Flange**



**1140x (CUT-TO-FIT) Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Convex
- ▶ Adhesive Border (with tape)

**HCPCS: A4407**




Stock No	Flange Size	For Stomas Up To
11402	 1-3/4" (44 mm)	1" (25 mm)
11403	 2-1/4" (57 mm)	1-1/2" (38 mm)
11404	 2-3/4" (70 mm)	2" (51 mm)



**1150x (PRE-SIZED) Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Convex
- ▶ Adhesive Border (with tape)

**HCPCS: A4407**

Stock No	Flange Size	Barrier Opening
11503	 1-3/4" (44 mm)	7/8" (22 mm)
11504	 1-3/4" (44 mm)	1" (25 mm)
11505	 2-1/4" (57 mm)	1-1/8" (29 mm)
11506	 2-1/4" (57 mm)	1-1/4" (32 mm)






**1530x (CUT-TO-FIT) Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Convex
- ▶ Full Hydrocolloid Barrier (without tape)

**HCPCS for 15302, 15303: A4407**

**HCPCS for 15304: A4408**

Stock No	Flange Size	For Stomas Up To
15302	 1-3/4" (44 mm)	1" (25 mm)
15303	 2-1/4" (57 mm)	1-1/2" (38 mm)
15304	 2-3/4" (70 mm)	2" (51 mm)

## CeraPlus™ Convex Skin Barriers with Integrated Floating Flange (continued)



**15310x** (PRE-SIZED) **Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Convex
- ▶ Full Hydrocolloid Barrier (without tape)



**HCPCS: A4407**

Stock No	Flange Size	Barrier Opening
153103	1-3/4" (44 mm)	7/8" (22 mm)
153104	1-3/4" (44 mm)	1" (25 mm)
153105	2-1/4" (57 mm)	1-1/8" (29 mm)
153106	2-1/4" (57 mm)	1-1/4" (32 mm)

## CeraPlus™ Flat Skin Barriers with Integrated Floating Flange



**1120x** (CUT-TO-FIT) **Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat
- ▶ Adhesive Border (with tape)

**HCPCS: A4409**

Stock No	Flange Size	For Stomas Up To
11202	1-3/4" (44 mm)	1-1/4" (32 mm)
11203	2-1/4" (57 mm)	1-3/4" (44 mm)
11204	2-3/4" (70 mm)	2-1/4" (57 mm)



**11201x** (PRE-SIZED) **Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat
- ▶ Adhesive Border (with tape)

**HCPCS: A4409**

Stock No	Flange Size	Barrier Opening
112013	1-3/4" (44 mm)	7/8" (22 mm)
112014	1-3/4" (44 mm)	1" (25 mm)
112015	1-3/4" (44 mm)	1-1/8" (29 mm)
112016	1-3/4" (44 mm)	1-1/4" (32 mm)
112017	2-1/4" (57 mm)	1-3/8" (35 mm)
112018	2-1/4" (57 mm)	1-1/2" (38 mm)



**1510x** (CUT-TO-FIT) **Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat
- ▶ Full Hydrocolloid Barrier (without tape)

**HCPCS for 15102, 15103: A4409**

**HCPCS for 15104: A4410**

Stock No	Flange Size	For Stomas Up To
15102	1-3/4" (44 mm)	1-1/4" (32 mm)
15103	2-1/4" (57 mm)	1-3/4" (44 mm)
15104	2-3/4" (70 mm)	2-1/4" (57 mm)



**15111x** (PRE-SIZED) **Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat
- ▶ Full Hydrocolloid Barrier (without tape)

**HCPCS: A4410**

Stock No	Flange Size	Barrier Opening
151113	1-3/4" (44 mm)	7/8" (22 mm)
151114	1-3/4" (44 mm)	1" (25 mm)
151115	1-3/4" (44 mm)	1-1/8" (29 mm)
151116	1-3/4" (44 mm)	1-1/4" (32 mm)
151117	2-1/4" (57 mm)	1-3/8" (35 mm)
151118	2-1/4" (57 mm)	1-1/2" (38 mm)

# Two-Piece Ostomy Pouches

## Available Features:

- Waterproof, integrated AF300™ filter is designed to allow gas and air to escape, but keeps odor and output in the pouch. Liquid cannot penetrate the filter, therefore not requiring a protective sticker
- Choice of secure pouch closure options – Lock 'n Roll™ microseal closure or clamp
- Numerous pouch size options to meet individual preferences
- Belt tabs for extra sense of security



## Drainable Pouches



### 1818x/1819x

Box of 10

- ▶ Drainable Pouch with *Lock 'n Roll*™ Microseal Closure
- ▶ Integrated *AF300*™ Filter
- ▶ Belt Tabs

#### HCPCS: A4425

Transparent Stock No	Beige Stock No	Flange Size	Pouch Length
18192	18182	1-3/4" (44 mm)	12" (30 cm)
18193	18183	2-1/4" (57 mm)	12" (30 cm)
18194	18184	2-3/4" (70 mm)	12" (30 cm)



### 1828x/1829x

Box of 20

- ▶ Drainable Mini-Pouch with *Lock 'n Roll*™ Microseal Closure
- ▶ Integrated *AF300*™ Filter
- ▶ Belt Tabs

#### HCPCS: A4425

Transparent Stock No	Beige Stock No	Flange Size	Pouch Length
18292	18282	1-3/4" (44 mm)	7" (18 cm)
18293	18283	2-1/4" (57 mm)	7" (18 cm)
18294	18284	2-3/4" (70 mm)	7" (18 cm)



### 1811x/1813x/1800x

Box of 10

- ▶ Drainable Pouch with *Lock 'n Roll*™ Microseal Closure
- ▶ Without Integrated Filter
- ▶ Belt Tabs

#### HCPCS: A5063

Transparent* Stock No	Beige Stock No	Flange Size	Pouch Length
18132	18112	1-3/4" (44 mm)	12" (30 cm)
18133	18113	2-1/4" (57 mm)	12" (30 cm)
18134	18114	2-3/4" (70 mm)	12" (30 cm)
18006	—	4" (102 mm)	12" (30 cm)



### 1814x/1816x

Box of 10  
Includes 1 Clamp

- ▶ Drainable Pouch with Clamp Closure
- ▶ Integrated *AF300*™ Filter
- ▶ Belt Tabs

#### HCPCS: A4425

Transparent Stock No	Beige Stock No	Flange Size	Pouch Length
18162	18142	1-3/4" (44 mm)	12" (30 cm)
18163	18143	2-1/4" (57 mm)	12" (30 cm)
18164	18144	2-3/4" (70 mm)	12" (30 cm)

## Drainable Pouches (continued)



### 1810x/1812x/1817x

Box of 10  
Includes 1 Clamp

- ▶ Drainable Pouch with Clamp Closure
- ▶ With Integrated *AF300*™ Filter
- ▶ Belt Tabs

#### HCPCS: A5063

Transparent Stock No	Beige Stock No	Flange Size	Pouch Length
18102	18122	1-3/4" (44 mm)	12" (30 cm)
18103	18123	2-1/4" (57 mm)	12" (30 cm)
18104	18124	2-3/4" (70 mm)	12" (30 cm)
18176	—	4" (102 mm)	12" (30 cm)



### 1820x

Box of 10  
Includes 1 Clamp

- ▶ Drainable Mini Pouch with Clamp Closure
- ▶ Without Integrated Filter
- ▶ Belt Tabs

#### HCPCS: A5063

Beige Stock No	Flange Size	Pouch Length
18202	1-3/4" (44 mm)	9" (23 cm)
18203	2-1/4" (57 mm)	9" (23 cm)
18204	2-3/4" (70 mm)	9" (23 cm)

## Closed Pouches



### 1836x/1837x

Box of 60

- ▶ Closed Pouch
- ▶ Integrated *AF300*™ Filter
- ▶ Belt Tabs

#### HCPCS: A4419

Transparent Stock No	Beige Stock No	Flange Size	Pouch Length
18362	18372	1-3/4" (44 mm)	9" (23 cm)
18363	18373	2-1/4" (57 mm)	9" (23 cm)
18364	18374	2-3/4" (70 mm)	9" (23 cm)



### 1838x/1839x

Box of 60

- ▶ Closed Mini-Pouch
- ▶ Integrated *AF300*™ Filter
- ▶ Belt Tabs

#### HCPCS: A4419

Transparent Stock No	Beige Stock No	Flange Size	Pouch Length
18382	18392	1-3/4" (44 mm)	7" (18 cm)
18383	18393	2-1/4" (57 mm)	7" (18 cm)
18384	18394	2-3/4" (70 mm)	7" (18 cm)

## Closed Pouches (continued)



**1835x**  
**1875x**

- ▶ Closed Mini-Pouch
- ▶ Without Integrated Filter
- ▶ Belt Tabs

**Box of 30**  
**Box of 60**

### HCPCS: A5054

Beige Stock No	Beige Stock No	Flange Size	Pouch Length
18352	18752	1-3/4" (44 mm)	7" (18 cm)
18353	18753	2-1/4" (57 mm)	7" (18 cm)
18354	—	2-3/4" (70 mm)	7" (18 cm)

## Urostomy Pouches



**1890x/1891x**

- ▶ Urostomy Pouch
- ▶ Soft, Streamline Tap
- ▶ Drain Valve Indicator
- ▶ Belt Tabs

**Box of 10**  
**Includes 2 adapters**

### HCPCS: A4432

Transparent Stock No	Beige Stock No	Flange Size	Pouch Length
18902	18912	1-3/4" (44 mm)	9" (23 cm)
18903	18913	2-1/4" (57 mm)	9" (23 cm)
18904	18914	2-3/4" (70 mm)	9" (23 cm)

## Non Sterile Kits — Drainable



**1960x**

**Box of 5**

Colostomy/Ileostomy Kit includes:

- ▶ (1) 1120x *CeraPlus*<sup>™</sup> (extended wear) Skin Barrier
- ▶ Flat Skin Barrier with Adhesive Border (with tape)
- ▶ (1) 1800x (*Lock 'n Roll*<sup>™</sup> Microseal Closure) Ultra-Clear Pouch Film
- ▶ Single-Use Package

### HCPCS: A5063 + A4414

Ultra-Clear Stock No	Flange Size	For Stomas Up To
19602	1-3/4" (44 mm)	1-1/4" (32 mm)
19603	2-1/4" (57 mm)	1-3/4" (44 mm)
19604	2-3/4" (70 mm)	2-1/4" (57 mm)



**1980x**

**Box of 5**

Colostomy/Ileostomy Kit includes:

- ▶ (1) 1120x *CeraPlus*<sup>™</sup> (extended wear) Skin Barrier
- ▶ Flat Skin Barrier with Adhesive Border (with tape)
- ▶ (1) 1817x (Clamp Closure) Ultra-Clear Pouch Film
- ▶ Single-Use Package
- ▶ Clamp Included

### HCPCS: A5063 + A4409

Ultra-Clear Stock No	Flange Size	For Stomas Up To
19802	1-3/4" (44 mm)	1-1/4" (32 mm)
19803	2-1/4" (57 mm)	1-3/4" (44 mm)
19804	2-3/4" (70 mm)	2-1/4" (57 mm)

## Non Sterile Kits — Urostomy



**1990x**

**Box of 5**

Urostomy Kit includes:

- ▶ (1) 1120x *CeraPlus*<sup>™</sup> (extended wear) Skin Barrier
- ▶ Flat Skin Barrier with Adhesive Border (with tape)
- ▶ (1) 1892x (Ultra-Clear Pouch Film)
- ▶ (1) 7331 Urostomy Drain Tube Adapter
- ▶ Single-Use Package

### HCPCS: A4432 + A4409

Ultra-Clear Stock No	Flange Size	For Stomas Up To
19902	1-3/4" (44 mm)	1-1/4" (32 mm)
19903	2-1/4" (57 mm)	1-3/4" (44 mm)
19904	2-3/4" (70 mm)	2-1/4" (57 mm)





## CeraPlus™ Soft Convex Skin Barrier - An Ideal Choice, Right from the Start

Choose a proactive path to healthy peristomal skin with CeraPlus™ Soft Convex Skin Barriers

### Meeting individual needs

CeraPlus™ Soft Convex Skin Barrier is available in a range of sizes in one- and two-piece to suit individual needs



**Comfortable and secure fit**  
Designed to conform to individual and changing body shapes for comfort and protection against leakage

### Built in belt tabs

The use of an ostomy belt can enhance the effect of the convexity by providing more barrier depth

### Ceramide-infused skin barrier

Ceramides serve as the skin's natural protection against water loss



### Integrated soft convex insert

Provides a secure and comfortable fit to support the peristomal skin while applying gentle pressure to help the stoma protrude and flatten creases and folds



Available in both one- and two-piece options.

“Convexity should be considered in the immediate post operative period to ensure a consistent predictable seal and reduce the risk of leakage.”

Convexity Consensus Statement 3  
Colwell J. et al. Use of a Convex Pouching System in the Post Operative Period. A National Consensus. JWOCN 2022;49(3):240-246

## CeraPlus™ One-Piece Pouching System

### Drainable Pouches with CeraPlus™ Soft Convex Skin Barriers



**895xx/896xx (CUT-TO-FIT)**

**Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Soft Convex
- ▶ Adhesive Border (with tape)
- ▶ Drainable Pouch with Lock 'n Roll™ Microseal Closure
- ▶ Integrated AF300™ Filter
- ▶ Belt Tabs

**HCPCS: A5057**

Beige w/ Viewing Option Stock No	Ultra-Clear Stock No	For Stomas Up To	Pouch Length
8954	8964	1" (25 mm)	12" (30 cm)
8958	8968	1-1/2" (38 mm)	12" (30 cm)
89511	89611	2-1/8" (55 mm)	12" (30 cm)



**89518/8195111/8195411 (CUT-TO-FIT)**

**Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Soft Convex
- ▶ Full Hydrocolloid Barrier (without tape)
- ▶ Drainable Pouch with Lock 'n Roll™ Microseal Closure
- ▶ Integrated AF300™ Filter
- ▶ Belt Tabs

**HCPCS: A5057**

Beige w/ Viewing Option Stock No	Ultra-Clear Stock No	For Stomas Up To	Pouch Length
89518	—	1-1/2" (38 mm)	12" (30 cm)
8195111	8195411	2-1/8" (55 mm)	12" (30 cm)



### DID YOU KNOW...

The CeraPlus™ Skin Barrier is infused with ceramide



- Ceramide is a natural component of the human skin
- Ceramide links the cells of the stratum corneum together to form a waterproof, protective barrier
- Ceramide helps prevent water loss that can lead to skin damage and dryness

Drainable Pouches with CeraPlus™ Soft Convex Skin Barriers (continued)



**896x** (PRE-SIZED) **Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Soft Convex
- ▶ Adhesive Border (with tape)
- ▶ Drainable Pouch with Lock 'n Roll™ Microseal Closure
- ▶ Integrated AF300™ Filter
- ▶ Belt Tabs

**HCPCS: A5057**

Beige w/Viewing Option Stock No	Barrier Opening	Pouch Length
8960	3/4" (19 mm)	12" (30 cm)
8961	7/8" (22 mm)	12" (30 cm)
8962	1" (25 mm)	12" (30 cm)
8963	1-1/8" (29 mm)	12" (30 cm)

Drainable Pouches with CeraPlus™ Convex Skin Barriers (continued)



**899x** (PRE-SIZED) **Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Convex
- ▶ Adhesive Border (with tape)
- ▶ Drainable Pouch with Lock 'n Roll™ Microseal Closure
- ▶ Integrated AF300™ Filter
- ▶ Belt Tabs

**HCPCS: A5057**

Beige w/Viewing Option Stock No	Barrier Opening	Pouch Length
8990	3/4" (19 mm)	12" (30 cm)
8991	7/8" (22 mm)	12" (30 cm)
8992	1" (25 mm)	12" (30 cm)
8993	1-1/8" (29 mm)	12" (30 cm)

Drainable Pouches with CeraPlus™ Convex Skin Barriers



**891xx** (CUT-TO-FIT) **Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Convex
- ▶ Adhesive Border (with tape)
- ▶ Drainable Pouch with Lock 'n Roll™ Microseal Closure
- ▶ Integrated AF300™ Filter
- ▶ Belt Tabs

**HCPCS: A5057**

Beige w/Viewing Option Stock No	For Stomas Up To	Pouch Length
8914	1" (25 mm)	12" (30 cm)
8918	1-1/2" (38 mm)	12" (30 cm)
89111	2" (51 mm)	12" (30 cm)

Drainable Pouches with CeraPlus™ Flat Skin Barriers



**8901/89321** (CUT-TO-FIT) **Box of 10**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat
- ▶ Adhesive Border (with tape)
- ▶ Drainable Pouch with Lock 'n Roll™ Microseal Closure
- ▶ Integrated AF300™ Filter

**HCPCS: A5056**

Beige w/Viewing Option Stock No	Transparent Stock No	For Stomas Up To	Pouch Length
8901	89321	2-1/2" (64 mm)	12" (30 cm)



**898xx** (CUT-TO-FIT) **Box of 5**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Convex
- ▶ Adhesive Border (with tape)
- ▶ Drainable Pouch with Lock 'n Roll™ Microseal Closure
- ▶ Without Integrated Filter
- ▶ Belt Tabs

**HCPCS: A4390**

Ultra-Clear Stock No	For Stomas Up To	Pouch Length
8984	1" (25 mm)	12" (30 cm)
8988	1-1/2" (38 mm)	12" (30 cm)
89811	2" (51 mm)	12" (30 cm)



**89311** (CUT-TO-FIT) **Box of 10**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat
- ▶ Full Hydrocolloid Barrier (without tape)
- ▶ Drainable Pouch with Lock 'n Roll™ Microseal Closure
- ▶ Integrated AF300™ Filter

**HCPCS: A5056**

Ultra-Clear Stock No	For Stomas Up To	Pouch Length
89311	2-1/2" (64 mm)	12" (30 cm)

## Drainable Pouches with CeraPlus™ Flat Skin Barriers (continued)



**8931** (CUT-TO-FIT) **Box of 10**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat
- ▶ Adhesive Border (with tape)
- ▶ Drainable Pouch with *Lock 'n Roll™* Microseal Closure
- ▶ Without Integrated Filter

**HCPCS: A4388**

Ultra-Clear Stock No	For Stomas Up To	Pouch Length
8931	2-1/2" (64 mm)	12" (30 cm)



**89301** (CUT-TO-FIT) **Box of 10**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat
- ▶ Full Hydrocolloid Barrier (without tape)
- ▶ Drainable Pouch with *Lock 'n Roll™* Microseal Closure
- ▶ Without Integrated Filter

**HCPCS: A5056**

Ultra-Clear Stock No	For Stomas Up To	Pouch Length
89301	2-1/2" (64 mm)	12" (30 cm)



**89xx** (PRE-SIZED) **Box of 10**

- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat
- ▶ Adhesive Border (with tape)
- ▶ Drainable Pouch with *Lock 'n Roll™* Microseal Closure
- ▶ Integrated *AF300™* Filter

**HCPCS: A5056**

Beige w/Viewing Option Stock No	Barrier Opening	Pouch Length
8925	1" (25 mm)	12" (30 cm)
8930	1-3/16" (30 mm)	12" (30 cm)
8935	1-3/8" (35 mm)	12" (30 cm)
8940	1-9/16" (40 mm)	12" (30 cm)

## Urostomy Pouches with CeraPlus™ Soft Convex Skin Barriers



SOFT CONVEX

**8413xx** (CUT-TO-FIT) **Box of 5**  
**Includes 1 adapter**

- ▶ Urostomy Pouch
- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Soft Convex
- ▶ Adhesive Border (with tape)
- ▶ Soft Streamline Tap
- ▶ Multi-Chamber Design
- ▶ Drain Valve Indicator
- ▶ Belt Tabs

**HCPCS: A4430**

Ultra-Clear Stock No	For Stomas Up To	Pouch Length
84134	1" (25 mm)	9" (23 cm)
84138	1-1/2" (38 mm)	9" (23 cm)
841311	2-1/8" (55 mm)	9" (23 cm)



SOFT CONVEX

**841x** (PRE-SIZED) **Box of 5**  
**Includes 1 adapter**

- ▶ Urostomy Pouch
- ▶ CeraPlus™ (extended wear) Skin Barrier
- ▶ Soft Convex
- ▶ Adhesive Border (with tape)
- ▶ Soft Streamline Tap
- ▶ Multi-Chamber Design
- ▶ Drain Valve Indicator
- ▶ Belt Tabs

**HCPCS: A4430**

Beige w/Viewing Option Stock No	Barrier Opening	Pouch Length
8412	3/4" (19 mm)	9" (23 cm)
8413	7/8" (22 mm)	9" (23 cm)
8414	1" (25 mm)	9" (23 cm)
8415	1-1/8" (29 mm)	9" (23 cm)





## Follow the evidence to help ensure healthy skin

After surgery: Peristomal skin starts healthy



Top layers of the epidermis are removed with every barrier change. It's unavoidable, and the damage isn't always visible.




Epidermal stripping is linked to increased levels of transepidermal water loss (TEWL)<sup>12</sup>

### References:

1. McNichol L, Lund C, Rosen T, Gray M. Medical adhesives and patient safety: state of the science: consensus statements for the assessment, prevention, and treatment of adhesive-related skin injuries. *J Wound Ostomy Continence Nurs.* 2013;40(4):365-380.
2. Grove G, Houser T, Sibbald G, Salvadalena G. Measuring epidermal effects of ostomy skin barriers. *Skin Res Technol.* 2019;25:179-186. <https://doi.org/10.1111/srt.12630>
3. Nybaek H, Lophagen S, Karlsmark T, Bang Knudsen D, Jemec GBE. Stratum corneum integrity as a predictor for peristomal skin problems in ostomates. *Br J Dermatol.* 2010;162:357-361.


## Healthy Skin: Keeps moisture in and irritants out

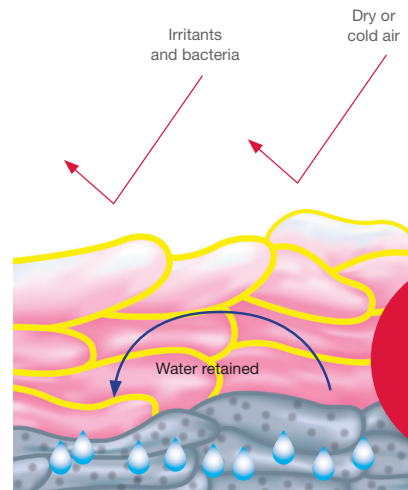

Ceramides are abundant in healthy skin. They form a waterproof barrier to help keep moisture in and bacteria out.

— = Ceramides

## Damaged Skin: Allows moisture out and irritants in

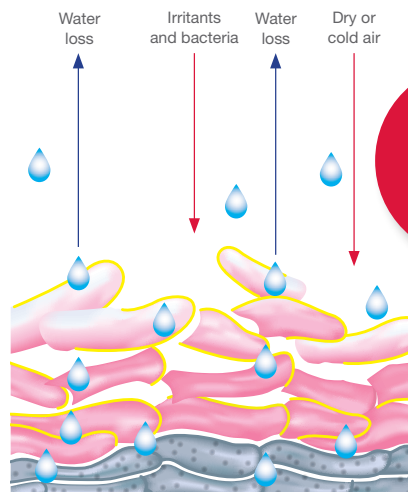

Damaged skin contains fewer ceramides.

💧 = Water droplet illustrates moisture loss from skin.

## The Results are Clear: CeraPlus™ Skin Barrier Helps Ensure Healthy Skin From Day One

At Hollister Ostomy Care, we believe peristomal skin deserves better. Clinical evidence continue to empower everything we do.

The ADVOCATE clinical trial looked at the difference in ostomy skin barriers and their impact on costs of stoma care and occurrence of peristomal skin complications (PSCs).

- 14% relative reduction in stoma-related cost of care with the CeraPlus™ Skin Barrier

### See the Evidence That Can Empower You.



CeraPlus™ Skin Barriers help support healthy skin and help minimize ostomy burden right from the start.<sup>1</sup>



CeraPlus™ Skin Barriers significantly reduce the stoma-related cost of care.<sup>1</sup>



Significantly more participants were "very satisfied" with overall barrier performance and itching prevention.<sup>1</sup>



15% fewer PSCs with the CeraPlus™ Skin Barrier.



	CeraPlus™ Skin Barrier	Coloplast SenSura Mio
Observed to be better at helping maintain intact skin. <sup>2,3</sup>		
Observed to be better at managing TEWL. <sup>2,3</sup>		
Observed to be better at protecting skin against damage from repeated barrier removal. <sup>2,3</sup>		

**Start using CeraPlus™ Skin Barrier from day one to help ensure the healthy skin your patients deserve.**

Contact your Hollister sales representative today for details.



References:

1. Colwell J, Pittman J, Salvadelena G. A randomized controlled trial detmrining variances in ostomy skin conditions and the economic impact (ADVOCATE). Journal Wound Ostomy Continence Nurs. 2018;45(1):37-42
2. Nichols T, Houser T, Grove G. Comparing the skin stripping effects of three ostomy skin barriers infused with ceramid, honey or aloe. Journal of Stomal Therapy Australia 2019; 39(2); 14-18
3. Grove G, Houser T, Sibbald G, Salvadelena G. Measuring epidermal effects of ostomy skin barriers. Skin Res Technol. 2019;25:179-186. <https://doi.org/10.1111/srt/12630>



Two-Piece  
Pouching System

# Two-Piece Pouching System

## Why choose a two-piece pouching system?

With a two-piece pouching system, the skin barrier and pouch are separate and connect together with a plastic ring, or flange. Here are a few reasons why you might want to use a two-piece pouching system for your patient:

- May want to change the pouch more often than the skin barrier
- May want to alternate between different pouch styles. For example, a smaller pouch may be worn during exercise or in the shower, and a larger pouch at night
- May find it easier to center the skin barrier over the stoma using a two-piece system





### Integrated Floating Flange

- Audible clicks signal a secure attachment
- Discreet low profile attachment of pouch to barrier



### CeraPlus™ Skin Barrier

CeraPlus™ Skin Barrier\* is infused with ceramide, the skin's naturally occurring protection against dryness. It is designed to maintain adhesive properties, and features a ceramide-infused formulation to help:

- Protect the skin's natural moisture barrier
- Maintain good peristomal skin health
- Decrease transepidermal water loss (TEWL) from damaged or eroded skin

*\*CeraPlus Skin Barriers contain the Remois technology of Alcare Co., Ltd.*

### Two-Piece Pouches



- Waterproof, integrated AF300™ filter is designed to allow gas and air to escape, but keeps odor and output in the pouch. Liquid cannot penetrate the filter, therefore not requiring a protective sticker
- Choice of secure pouch closure options – clamp or integrated closure
- Numerous pouch size options to meet individual preferences

## Two-Piece Pouching System

For corresponding pouches, please see pages 18-20.

### Flextend™ Convex Skin Barriers



Box of 5

(CUT-TO-FIT)

Stock No	Flange Size	For Stomas Up To	Tape/ No Tape	CeraPlus™ Product		Page
				HCPCS	Alternative	
14802	1-3/4" (44mm)	1" (25mm)	TAPE	A4407	11402	15
14803	2-1/4" (57mm)	1-1/2" (38mm)	TAPE	A4407	11403	15
14804	2-3/4" (70mm)	2" (51mm)	TAPE	A4407	11404	15
15802	1-3/4" (44mm)	1" (25mm)	NO TAPE	A4407	15302	15
15803	2-1/4" (57mm)	1-1/2" (38mm)	NO TAPE	A4407	15303	15
15804	2-3/4" (70mm)	2" (51mm)	NO TAPE	A4408	15304	15

(PRE-SIZED)

Stock No	Flange Size	Barrier Opening	Tape/ No Tape	CeraPlus™ Product		Page
				HCPCS	Alternative	
13901	1-3/4" (44mm)	5/8" (16mm)	TAPE	A4407	—	—
13902	1-3/4" (44mm)	3/4" (19mm)	TAPE	A4407	—	—
13903	1-3/4" (44mm)	7/8" (22mm)	TAPE	A4407	11503	15
13904	1-3/4" (44mm)	1" (25mm)	TAPE	A4407	11504	15
13905	2-1/4" (57mm)	1-1/8" (29mm)	TAPE	A4407	11505	15
13906	2-1/4" (57mm)	1-1/4" (32mm)	TAPE	A4407	11506	15
13907	2-1/4" (57mm)	1-3/8" (35mm)	TAPE	A4407	—	—
13908	2-1/4" (57mm)	1-1/2" (38mm)	TAPE	A4407	—	—
13909	2-3/4" (70mm)	1-5/8" (41mm)	TAPE	A4407	—	—
13910	2-3/4" (70mm)	1-3/4" (44mm)	TAPE	A4407	—	—
13911	2-3/4" (70mm)	2" (51mm)	TAPE	A4407	—	—
15903	1-3/4" (44mm)	7/8" (22mm)	NO TAPE	A4407	153103	15
15904	1-3/4" (44mm)	1" (25mm)	NO TAPE	A4407	153104	15
15905	2-1/4" (57mm)	1-1/8" (29mm)	NO TAPE	A4407	153105	15
15906	2-1/4" (57mm)	1-1/4" (32mm)	NO TAPE	A4407	153106	15

CeraPlus™ Product Alternative listed may have slightly different features, please see corresponding page for details.

For corresponding pouches, please see pages 18-20.

### Flextend™ Flat Skin Barriers

Box of 5

(CUT-TO-FIT)

Stock No	Flange Size	For Stomas Up To	Tape/ No Tape	CeraPlus™ Product		Page
				HCPCS	Alternative	
14602	1-3/4" (44mm)	1-1/4" (32mm)	TAPE	A4409	11202	16
14603	2-1/4" (57mm)	1-3/4" (44mm)	TAPE	A4409	11203	16
14604	2-3/4" (70mm)	2-1/4" (57mm)	TAPE	A4409	11204	16
14606	4" (102mm)	3-1/2" (89mm)	TAPE	A4410	—	—
15602	1-3/4" (44mm)	1-1/4" (32mm)	NO TAPE	A4409	15102	16
15603	2-1/4" (57mm)	1-3/4" (44mm)	NO TAPE	A4409	15103	16
15604	2-3/4" (70mm)	2-1/4" (57mm)	NO TAPE	A4410	15104	16

(PRE-SIZED)

Stock No	Flange Size	Barrier Opening	Tape/ No Tape	CeraPlus™ Product		Page
				HCPCS	Alternative	
14701	1-3/4" (44mm)	5/8" (16mm)	TAPE	A4409	—	—
14702	1-3/4" (44mm)	3/4" (19mm)	TAPE	A4409	—	—
14703	1-3/4" (44mm)	7/8" (22mm)	TAPE	A4409	112013	16
14704	1-3/4" (44mm)	1" (25mm)	TAPE	A4409	112014	16
14705	1-3/4" (44mm)	1-1/8" (29mm)	TAPE	A4409	112015	16
14706	1-3/4" (44mm)	1-1/4" (32mm)	TAPE	A4409	112016	16
14707	2-1/4" (57mm)	1-3/8" (35mm)	TAPE	A4409	112017	16
14708	2-1/4" (57mm)	1-1/2" (38mm)	TAPE	A4409	112018	16
14709	2-1/4" (57mm)	1-3/4" (44mm)	TAPE	A4409	—	—
16101	1-3/4" (44mm)	5/8" (16mm)	NO TAPE	A4409	—	—
16102	1-3/4" (44mm)	3/4" (19mm)	NO TAPE	A4409	—	—
16103	1-3/4" (44mm)	7/8" (22mm)	NO TAPE	A4409	151113	16
16104	1-3/4" (44mm)	1" (25mm)	NO TAPE	A4409	151114	16
16105	1-3/4" (44mm)	1-1/8" (29mm)	NO TAPE	A4409	151115	16
16106	1-3/4" (44mm)	1-1/4" (32mm)	NO TAPE	A4409	151116	16
16107	2-1/4" (57mm)	1-3/8" (35mm)	NO TAPE	A4409	151117	16
16108	2-1/4" (57mm)	1-1/2" (38mm)	NO TAPE	A4409	151118	16

CeraPlus™ Product Alternative listed may have slightly different features, please see corresponding page for details.







## Two-Piece Pouching System

For corresponding pouches, please see pages 18-20.









### FlexWear™ Convex Skin Barriers

 CONVEX **Box of 5**

(CUT-TO-FIT)

Stock No	Flange Size	For Stomas Up To	Tape/ No Tape	HCPCS	CeraPlus™ Product Alternative	Page
13402	 1-3/4" (44mm)	1" (25mm)	TAPE	A4373	11402	15
13403	 2-1/4" (57mm)	1-1/2" (38mm)	TAPE	A4373	11403	15
13404	 2-3/4" (70mm)	2" (51mm)	TAPE	A4373	11404	15
15402	 1-3/4" (44mm)	1" (25mm)	NO TAPE	A4373	15302	15
15403	 2-1/4" (57mm)	1-1/2" (38mm)	NO TAPE	A4373	15303	15
15404	 2-3/4" (70mm)	2" (51mm)	NO TAPE	A4373	15304	15







(PRE-SIZED)

Stock No	Flange Size	Barrier Opening	Tape/ No Tape	HCPCS	CeraPlus™ Product Alternative	Page
13503	 1-3/4" (44mm)	7/8" (22mm)	TAPE	A4373	11503	15
13504	 1-3/4" (44mm)	1" (25mm)	TAPE	A4373	11504	15
13505	 2-1/4" (57mm)	1-1/8" (29mm)	TAPE	A4373	11505	15
13506	 2-1/4" (57mm)	1-1/4" (32mm)	TAPE	A4373	11506	15
15503	 1-3/4" (44mm)	7/8" (22mm)	NO TAPE	A4373	153103	16
15504	 1-3/4" (44mm)	1" (25mm)	NO TAPE	A4373	153104	16
15505	 2-1/4" (57mm)	1-1/8" (29mm)	NO TAPE	A4373	153105	16
15506	 2-1/4" (57mm)	1-1/4" (32mm)	NO TAPE	A4373	153106	16

















### FlexWear™ Flat Skin Barriers

**Box of 5**

(CUT-TO-FIT)

Stock No	Flange Size	For Stomas Up To	Tape/ No Tape	HCPCS	CeraPlus™ Product Alternative	Page
14202	 1-3/4" (44mm)	1-1/4" (32mm)	TAPE	A4414	11202	16
14203	 2-1/4" (57mm)	1-3/4" (57mm)	TAPE	A4414	11203	16
14204	 2-3/4" (70mm)	2-1/4" (57mm)	TAPE	A4414	11204	16
15202	 1-3/4" (44mm)	1-1/4" (32mm)	NO TAPE	A4414	15102	16
15203	 2-1/4" (57mm)	1-3/4" (57mm)	NO TAPE	A4414	15103	16
15204	 2-3/4" (70mm)	2-1/4" (57mm)	NO TAPE	A4415	15104	16

(PRE-SIZED)

Stock No	Flange Size	Barrier Opening	Tape/ No Tape	HCPCS	CeraPlus™ Product Alternative	Page
14301	 1-3/4" (44mm)	5/8" (16mm)	TAPE	A4414	—	—
14302	 1-3/4" (44mm)	3/4" (19mm)	TAPE	A4414	—	—
14303	 1-3/4" (44mm)	7/8" (22mm)	TAPE	A4414	112013	16
14304	 1-3/4" (44mm)	1" (25mm)	TAPE	A4415	112014	16
14305	 1-3/4" (44mm)	1-1/8" (29mm)	TAPE	A4414	112015	16
14306	 1-3/4" (44mm)	1-1/4" (32mm)	TAPE	A4414	112016	16
14307	 2-1/4" (57mm)	1-3/8" (35mm)	TAPE	A4414	112017	16
14308	 2-1/4" (57mm)	1-1/2" (38mm)	TAPE	A4414	112018	16
14309	 2-1/4" (57mm)	1-3/4" (44mm)	TAPE	A4414	—	—
16402	 1-3/4" (44mm)	3/4" (19mm)	NO TAPE	A4414	—	—
16403	 1-3/4" (44mm)	7/8" (22mm)	NO TAPE	A4414	151113	16
16404	 1-3/4" (44mm)	1" (25mm)	NO TAPE	A4414	151114	16
16405	 1-3/4" (44mm)	1-1/8" (29mm)	NO TAPE	A4414	151115	16
16406	 1-3/4" (44mm)	1-1/4" (32mm)	NO TAPE	A4414	151116	16
16407	 2-1/4" (57mm)	1-3/8" (35mm)	NO TAPE	A4414	151117	16
16408	 2-1/4" (57mm)	1-1/2" (38mm)	NO TAPE	A4414	151118	16

*CeraPlus™ Product Alternative listed may have slightly different features, please see corresponding page for details.*



One-Piece  
Pouching System

One-Piece  
Pouching System

## Why choose a one-piece pouching system?

With a one-piece pouching system, the skin barrier and pouch are connected and cannot be separated. Here are a few reasons why you might want to use a one-piece pouching system:

- Need a pouching system that is flexible
- May prefer a pouching system that does not require putting pieces together
- Do not wish to change the pouch more often than the skin barrier

# One-Piece Pouching System

## Product Features Key

reference for products listed on pages 38-42



**UCLR:** Ultra-Clear

---



**BGVO:** Beige with viewing option

---



**BG:** Beige

---



**LNR:** Drainable pouch with  
*Lock 'N Roll'* (integrated) Closure

---



**CLAMP:** Drainable pouch with clamp closure

---



**URO:** Urostomy pouching system

---



**CLOSED:** Closed pouching system

## Flexlend™ Soft Convex Skin Barriers



All One-Piece Soft Convex Pouching Systems include Belt Tabs (CUT-TO-FIT)

- ▶ Adhesive Border (with tape)
- ▶ With Integrated AF300™ Filter\*

Stock No	Color	For Stomas Up To	Length	Closure Type	HPCPS	CeraPlus™ Product Alternative	Page
8574*	UCLR	1" (25mm)	12" (30cm)	LNR	A5057	8964	23
8674*	BGVO	1" (25mm)	12" (30cm)	LNR	A5057	8954	23
8578*	UCLR	1-1/2" (38mm)	12" (30cm)	LNR	A5057	8968	23
8678*	BGVO	1-1/2" (38mm)	12" (30cm)	LNR	A5057	8958	23
85711*	UCLR	2-1/8" (55mm)	12" (30cm)	LNR	A5057	89611	23
86711*	BGVO	2-1/8" (55mm)	12" (30cm)	LNR	A5057	89511	23

(PRE-SIZED)

- ▶ Adhesive Border (with tape)
- ▶ With Integrated AF300™ Filter\*

Stock No	Color	Barrier Opening	Length	Closure Type	HPCPS	CeraPlus™ Product Alternative	Page
8661*	BGVO	7/8" (22mm)	12" (30cm)	LNR	A5057	8961	24
8662*	BGVO	1" (25mm)	12" (30cm)	LNR	A5057	8962	24
8663*	BGVO	1-1/8" (29mm)	12" (30cm)	LNR	A5057	8963	24
8664*	BGVO	1-1/4" (32mm)	12" (30cm)	LNR	A5057	—	—
8665*	BGVO	1-1/2" (38mm)	12" (30cm)	LNR	A5057	—	—
8666*	BGVO	1-3/4" (44mm)	12" (30cm)	LNR	A5057	—	—

CeraPlus™ Product Alternative listed may have slightly different features, please see corresponding page for details.

## Flexlend™ Convex Skin Barriers



All One-Piece Convex Pouching Systems include Belt Tabs (CUT-TO-FIT)

- ▶ Adhesive Border (with tape)

Stock No	Color	For Stomas Up To	Length	Closure Type	HPCPS	CeraPlus™ Product Alternative	Page
8274*	UCLR	1" (25mm)	12" (30cm)	LNR	A5057	8914	24
8544*	BGVO	1" (25mm)	12" (30cm)	LNR	A5057	8914	24
8278*	UCLR	1-1/2" (38mm)	12" (30cm)	LNR	A5057	8918	24
8548*	BGVO	1-1/2" (38mm)	12" (30cm)	LNR	A5057	8918	24
82711*	UCLR	2" (51mm)	12" (30cm)	LNR	A5057	89111	24
8244	UCLR	1" (25mm)	12" (30cm)	LNR	A4390	8984	24
8248	UCLR	1-1/2" (38mm)	12" (30cm)	LNR	A4390	8988	24
82411	UCLR	2" (51mm)	12" (30cm)	LNR	A4390	89811	24
82624	UCLR	1" (25mm)	12" (30cm)	CLAMP	A4390	8984	24
82628	UCLR	1-1/2" (38mm)	12" (30cm)	CLAMP	A4390	8988	24
826211	UCLR	2" (51mm)	12" (30cm)	CLAMP	A4390	89811	24
84794	UCLR	1" (25mm)	9" (23cm)	URO	A4430	84404	28
84798	UCLR	1-1/2" (38mm)	9" (23cm)	URO	A4430	84408	28
8479111	UCLR	2" (51mm)	9" (23cm)	URO	A4430	844011	28
8439111	BG	2" (51mm)	9" (23cm)	URO	A4430	—	—

With Integrated AF300™ Filter\*

URO pouches come with 1 adapter/box

CLAMP pouches come with 1 clamp/box

CeraPlus™ Product Alternative listed may have slightly different features, please see corresponding page for details.

## Flexlend™ Convex Skin Barriers (continued)

All One-Piece Convex Pouching Systems include Belt Tabs (PRE-SIZED)

- ▶ Adhesive Border (with tape)

Stock No	Color	Barrier Opening	Length	Closure Type	HPCPS	CeraPlus™ Product Alternative	Page
8280	UCLR	3/4" (19mm)	12" (30cm)	LNR	A4390	8990	25
8290	BG	3/4" (19mm)	12" (30cm)	LNR	A4390	8990	25
8281	UCLR	7/8" (22mm)	12" (30cm)	LNR	A4390	8991	25
8291	BG	7/8" (22mm)	12" (30cm)	LNR	A4390	8991	25
8282	UCLR	1" (25mm)	12" (30cm)	LNR	A4390	8992	25
8292	BG	1" (25mm)	12" (30cm)	LNR	A4390	8992	25
8283	UCLR	1-1/8" (29mm)	12" (30cm)	LNR	A4390	8993	25
8284	UCLR	1-1/4" (32mm)	12" (30cm)	LNR	A4390	—	—
8294	BG	1-1/4" (32mm)	12" (30cm)	LNR	A4390	—	—
8285	UCLR	1-3/8" (35mm)	12" (30cm)	LNR	A4390	—	—
8295	BG	1-3/8" (35mm)	12" (30cm)	LNR	A4390	—	—
8286	UCLR	1-1/2" (38mm)	12" (30cm)	LNR	A4390	—	—
8288	UCLR	1-3/4" (44mm)	12" (30cm)	LNR	A4390	—	—
82610	UCLR	3/4" (19mm)	12" (30cm)	CLAMP	A4390	8990	25
82611	UCLR	7/8" (22mm)	12" (30cm)	CLAMP	A4390	8991	25
82612	UCLR	1" (25mm)	12" (30cm)	CLAMP	A4390	8992	25
82613	UCLR	1-1/8" (29mm)	12" (30cm)	CLAMP	A4390	8993	25
82614	UCLR	1-1/4" (32mm)	12" (30cm)	CLAMP	A4390	—	—
82615	UCLR	1-3/8" (35mm)	12" (30cm)	CLAMP	A4390	—	—
82616	UCLR	1-1/2" (38mm)	12" (30cm)	CLAMP	A4390	—	—
82617	UCLR	1-5/8" (41mm)	12" (30cm)	CLAMP	A4390	—	—
82619	UCLR	2" (51mm)	12" (30cm)	CLAMP	A4390	—	—
84890	UCLR	1/2" (13mm)	9" (23cm)	URO	A4430	—	—
84891	UCLR	5/8" (16mm)	9" (23cm)	URO	A4430	—	—
84991	BG	5/8" (16mm)	9" (23cm)	URO	A4430	—	—
84892	UCLR	3/4" (19mm)	9" (23cm)	URO	A4430	84292	28
84992	BG	3/4" (19mm)	9" (23cm)	URO	A4430	84292	28
84893	UCLR	7/8" (22mm)	9" (23cm)	URO	A4430	84293	28
84993	BG	7/8" (22mm)	9" (23cm)	URO	A4430	84293	28
84894	UCLR	1" (25mm)	9" (23cm)	URO	A4430	84294	28
84994	BG	1" (25mm)	9" (23cm)	URO	A4430	84294	28
84895	UCLR	1-1/8" (29mm)	9" (23cm)	URO	A4430	84295	28
84995	BG	1-1/8" (29mm)	9" (23cm)	URO	A4430	84295	28
84896	UCLR	1-1/4" (32mm)	9" (23cm)	URO	A4430	—	—
84996	BG	1-1/4" (32mm)	9" (23cm)	URO	A4430	—	—
84897	UCLR	1-3/8" (35mm)	9" (23cm)	URO	A4430	—	—
84997	BG	1-3/8" (35mm)	9" (23cm)	URO	A4430	—	—
84898	UCLR	1-1/2" (38mm)	9" (23cm)	URO	A4430	—	—
84998	BG	1-1/2" (38mm)	9" (23cm)	URO	A4430	—	—
84899	UCLR	1-3/4" (44mm)	9" (23cm)	URO	A4430	—	—
84999	BG	1-3/4" (44mm)	9" (23cm)	URO	A4430	—	—

URO pouches come with 1 adapter/box

CLAMP pouches come with 1 clamp/box

CeraPlus™ Product Alternative listed may have slightly different features, please see corresponding page for details.

## Flextend™ Flat Skin Barriers

Box of 10

(CUT-TO-FIT)

- ▶ Adhesive Border (with tape)

Stock No	Color	For Stomas Up To	Length	Closure Type	HCPCS	CeraPlus™ Product	
						Alternative	Page
8331*	UCLR	2-1/2" (64mm)	12" (30cm)	LNR	A5056	89321	25
8371*	BGVO	2-1/2" (64mm)	12" (30cm)	LNR	A5056	8901	25
8531	UCLR	2-1/2" (64mm)	12" (30cm)	LNR	A5056	8931	26
8551	BG	2-1/2" (64mm)	12" (30cm)	LNR	A5056	8931	26
8631	UCLR	2-1/2" (64mm)	12" (30cm)	CLAMP	A5056	8931	26
8641	UCLR	2-1/2" (64mm)	9" (23cm)	CLAMP	A5056	—	—
84590	UCLR	2-1/2" (64mm)	9" (23cm)	URO	A4428	—	—
84490	BG	2-1/2" (64mm)	9" (23cm)	URO	A4428	—	—

With Integrated AF300™ Filter\*

URO pouches come with 2 adapters/box

CLAMP pouches come with 1 clamp/box

(PRE-SIZED)

- ▶ Adhesive Border (with tape)

Stock No	Color	Barrier Opening	Length	Closure Type	HCPCS	CeraPlus™ Product	
						Alternative	Page
8552	BG	1" (25mm)	12" (30cm)	LNR	A4388	8925	26
8538	UCLR	1-1/4" (32mm)	12" (30cm)	LNR	A4388	—	—
8558	BG	1-1/4" (32mm)	12" (30cm)	LNR	A4388	—	—
8637	UCLR	3/4" (19mm)	12" (30cm)	CLAMP	A4388	—	—
8632	UCLR	1" (25mm)	12" (30cm)	CLAMP	A4388	8925	26
8638	UCLR	1-1/4" (32mm)	12" (30cm)	CLAMP	A4388	—	—
8633	UCLR	1-1/2" (38mm)	12" (30cm)	CLAMP	A4388	—	—
8639	UCLR	1-3/4" (44mm)	12" (30cm)	CLAMP	A4388	—	—
8634	UCLR	2" (51mm)	12" (30cm)	CLAMP	A4388	—	—
8635	UCLR	2-1/2" (64mm)	12" (30cm)	CLAMP	A4388	—	—
8647	BG	3/4" (19mm)	9" (23cm)	CLAMP	A4388	—	—
8642	BG	1" (25mm)	9" (23cm)	CLAMP	A4388	—	—
8648	BG	1-1/4" (32mm)	9" (23cm)	CLAMP	A4388	—	—
8643	BG	1-1/2" (38mm)	9" (23cm)	CLAMP	A4388	—	—
8649	BG	1-3/4" (44mm)	9" (23cm)	CLAMP	A4388	—	—
84497	BG	3/4" (19mm)	9" (23cm)	URO	A4428	—	—
84692	UCLR	1" (25mm)	9" (23cm)	URO	A4428	—	—
84698	UCLR	1-1/4" (32mm)	9" (23cm)	URO	A4428	—	—
84498	BG	1-1/4" (32mm)	9" (23cm)	URO	A4428	—	—
84693	UCLR	1-1/2" (38mm)	9" (23cm)	URO	A4428	—	—
84699	UCLR	1-3/4" (44mm)	9" (23cm)	URO	A4428	—	—
84499	BG	1-3/4" (44mm)	9" (23cm)	URO	A4428	—	—
84694	UCLR	2" (51mm)	9" (23cm)	URO	A4428	—	—
84494	BG	2" (51mm)	9" (23cm)	URO	A4428	—	—

URO pouches come with 2 adapters/box

CLAMP pouches come with 1 clamp/box

CeraPlus™ Product Alternative listed may have slightly different features, please see corresponding page for details.

## FlexWear™ Flat Skin Barriers

Box of 10

(CUT-TO-FIT)

- ▶ Adhesive Border (with tape)
- ▶ With Integrated AF300™ Filter

Stock No	Color	For Stomas Up To	Length	Closure Type	HCPCS	CeraPlus™ Product	
						Alternative	Page
8081*	UCLR	2-1/2" (64mm)	12" (30cm)	LNR	A4424	89321	25
8181*	BG	2-1/2" (64mm)	12" (30cm)	LNR	A4424	8901	25
8381*	BGVO	2-1/2" (64mm)	12" (30cm)	LNR	A4424	8901	25
8031	UCLR	2-1/2" (64mm)	12" (30cm)	LNR	A5061	8931	25

With Integrated AF300™ Filter\*

URO pouches come with 1 adapters/box

CeraPlus™ Product Alternative listed may have slightly different features, please see corresponding page for details.



# One-Piece Pouching System

## SoftFlex™ Flat Skin Barriers

(CUT-TO-FIT)

- ▶ Full Hydrocolloid Barrier (without tape)

883x/884x/885x: Box of 10

887x/888x: Box of 20

821x/823x/824x: Box of 30

Stock No	Color	For Stomas Up To	Length	Closure Type	HCPCS	CeraPlus™ Product Alternative	Page
88400*	UCLR	2-1/8" (55mm)	12" (30cm)	LNR	A4424	89311	25
88300*	BG	2-1/8" (55mm)	12" (30cm)	LNR	A4424	89311	25
88800*	UCLR	2-1/8" (55mm)	7" (18cm)	LNR	A4424	—	—
88700*	BG	2-1/8" (55mm)	7" (18cm)	LNR	A4424	—	—
88402*	UCLR	3" x 2.5" (75mm x 65mm)	12" (30cm)	LNR	A4424	—	—
88302*	BG	3" x 2.5" (75mm x 65mm)	12" (30cm)	LNR	A4424	—	—
88501	UCLR	2-1/2" (64mm)	12" (30cm)	CLAMP	A5061	89301	26
82400*	UCLR	2-1/8" (55mm)	9" (23cm)	CLOSED	A4416	—	—
82300*	BG	2-1/8" (55mm)	9" (23cm)	CLOSED	A4416	—	—
82402*	UCLR	3" x 2.5" (75mm x 65mm)	9" (23cm)	CLOSED	A4416	—	—
82302*	BG	3" x 2.5" (75mm x 65mm)	9" (23cm)	CLOSED	A4416	—	—
82500*	UCLR	2-1/8" (55mm)	7" (18cm)	CLOSED	A4416	—	—
82100*	BG	2-1/8" (55mm)	7" (18cm)	CLOSED	A4416	—	—

With Integrated AF300™ Filter\*

(PRE-CUT)

- ▶ Full Hydrocolloid Barrier (without tape)

Stock No	Color	Barrier Opening	Length	Closure Type	HCPCS	CeraPlus™ Product Alternative	Page
88425*	UCLR	1" (25mm)	12" (30cm)	LNR	A4424	—	—
88325*	BG	1" (25mm)	12" (30cm)	LNR	A4424	—	—
88725*	BG	1" (25mm)	7" (18cm)	LNR	A4424	—	—
88430*	UCLR	1-3/16" (30mm)	12" (30cm)	LNR	A4424	—	—
88330*	BG	1-3/16" (30mm)	12" (30cm)	LNR	A4424	—	—
88730*	BG	1-3/16" (30mm)	7" (18cm)	LNR	A4424	—	—
88435*	UCLR	1-3/8" (35mm)	12" (30cm)	LNR	A4424	—	—
88335*	BG	1-3/8" (35mm)	12" (30cm)	LNR	A4424	—	—
88735*	BG	1-3/8" (35mm)	7" (18cm)	LNR	A4424	—	—
88440*	UCLR	1-9/16" (40mm)	12" (30cm)	LNR	A4424	—	—
88340*	BG	1-9/16" (40mm)	12" (30cm)	LNR	A4424	—	—
88740*	BG	1-9/16" (40mm)	7" (18cm)	LNR	A4424	—	—
82325*	BG	1" (25mm)	9" (23cm)	CLOSED	A4416	—	—
82330*	BG	1-3/16" (30mm)	9" (23cm)	CLOSED	A4416	—	—
82335*	BG	1-3/8" (35mm)	9" (23cm)	CLOSED	A4416	—	—
82125*	BG	1" (25mm)	7" (18cm)	CLOSED	A4416	—	—
82130*	BG	1-3/16" (30mm)	7" (18cm)	CLOSED	A4416	—	—
82135*	BG	1-3/8" (35mm)	7" (18cm)	CLOSED	A4416	—	—

With Integrated AF300™ Filter\*

CeraPlus™ Product Alternative listed may have slightly different features, please see corresponding page for details.



## A Flexible Fit Made Secure

### CeraPlus™ Soft Convex Pouching System

Designed to stay secure, no matter what shape the body takes.

When it comes to finding the right solution for patients, fit comes first. Comfort and ease-of-use are also essential to help ensure positive outcomes—in the hospital and after patients leave your care.

CeraPlus™ Soft Convex Skin Pouching System delivers both convex support and flexibility. The skin barrier is designed to conform to the body, without compromising fit. Belt tabs and a tape-bordered barrier help provide an extra sense of security.

- Look for the soft convex symbol 





Hospital Products

Hospital Products



## Two-Piece High Output






### 1802x

Box of 10

- ▶ Drainable Pouch
  - ▶ Integrated AF300™ Filter
  - ▶ Wide Bore Tap (16.5 mm inner diameter)
  - ▶ Belt Tabs
  - ▶ No adapter needed to connect to drainage collector 5551.\*  
Does not fit a standard urinary bedside collector
- \* See page 62 for bedside drainage collector 5551

#### HCPCS: A4413

Ultra-Clear Stock No	Flange Size	Pouch Length
18022	 1-3/4" (44 mm)	12" (30 cm)
18023	 2-1/4" (57 mm)	12" (30 cm)
18024	 2-3/4" (70 mm)	12" (30 cm)







### 1801x

Box of 10

- ▶ Drainable Pouch
  - ▶ Without Integrated Filter
  - ▶ Soft Drain Tap
  - ▶ Belt Tabs
  - ▶ No adapter needed to connect to bedside collector 9839.
- \* See page 62 for bedside drainage collector 9839

#### HCPCS: A4412

Ultra-Clear Stock No	Flange Size	Pouch Length
18012	 1-3/4" (44 mm)	12" (30 cm)
18013	 2-1/4" (57 mm)	12" (30 cm)
18014	 2-3/4" (70 mm)	12" (30 cm)
18016	 4" (102mm)	12" (30 cm)

## One-Piece High Output

### Flexend™ Soft Convex Skin Barriers

Box of 10



(CUT-TO-FIT)

- ▶ Adhesive Border (with tape)
- ▶ With Integrated AF300® Filter
- ▶ Wide Bore Tap (16.5 mm inner diameter)
- ▶ Belt Tabs
- ▶ No adapter needed to connect to drainage collector 5551.\*  
Does not fit a standard urinary bedside collector

\* See page 62 for bedside drainage collector 5551

Stock No	Color	For Stomas Up To	Length	Tap Type	HCPCS
83025	UCLR	1" (25mm)	12" (30cm)	WIDE	A4435
83038	UCLR	1-1/2" (38mm)	12" (30cm)	WIDE	A4435
83055	UCLR	2-1/8" (55mm)	12" (30cm)	WIDE	A4435

### Flexend™ Flat Skin Barriers



(CUT-TO-FIT)

- ▶ Adhesive Border (with tape)
- ▶ With Integrated AF300® Filter
- ▶ Wide Bore Tap (16.5 mm inner diameter)
- ▶ Belt Tabs
- ▶ No adapter needed to connect to drainage collector 5551.\*  
Does not fit a standard urinary bedside collector

\* See page 62 for bedside drainage collector 5551

Stock No	Color	For Stomas Up To	Length	Tap Type	HCPCS
26064	UCLR	2-1/2" (64mm)	12" (30cm)	WIDE	A4435



(CUT-TO-FIT)

- ▶ Full Hydrocolloid Barrier (without tape)
- ▶ No Starter Hole
- ▶ Access Window
- ▶ No adapter needed to connect to bedside collector 9839.

\* See page 62 for bedside drainage collector 9839

Stock No	Color	For Stomas Up To	Length	Closure Type	HCPCS
80070	UCLR	2-3/4" (70mm)	12" (30cm)	SOFT	A4435
80110	UCLR	4-1/3" (110mm)	12" (30cm)	SOFT	A4435

## Wound Drainage Collectors



- ▶ Full Hydrocolloid barrier (without tape)
  - ▶ No Starter hole
  - ▶ Access Window
  - ▶ No adapter needed to connect to bedside collector 9839.
- \* See page 62 for bedside drainage collector 9839

Stock No	Type	Formulation	Description	Qty/Box	HCPCS
9773	Non Sterile	FlexWear	Small (for wounds up to 3")	10	A6154
9776	Non Sterile	FlexWear	Medium (for wounds up to 3-3/4")	10	A6154
9778	Non Sterile	FlexWear	Large (for wounds up to 4" x 8")	3	A6154
9772	Sterile	No Barrier	Small (for wounds up to 3")	10	A6154
9775	Sterile	No Barrier	Medium (for wounds up to 3-3/4")	10	A6154

## Hospital Kits

## Non Sterile Kits — Drainable



Ceramide

## 1960x

Box of 5

Colostomy/Ileostomy Kit includes:

- ▶ (1) 1120x CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat Skin Barrier with Adhesive Border (with tape)
- ▶ (1) 1800x (Lock 'n Roll™ Microseal Closure) Ultra-Clear Pouch Film
- ▶ Single-Use Package

HCPCS: A5063 + A4414

Ultra-Clear Stock No	Flange Size	For Stomas Up To
19602	1-3/4" (44 mm)	1-1/4" (32 mm)
19603	2-1/4" (57 mm)	1-3/4" (44 mm)
19604	2-3/4" (70 mm)	2-1/4" (57 mm)



Ceramide

## 1980x

Box of 5

Colostomy/Ileostomy Kit includes:

- ▶ (1) 1120x CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat Skin Barrier with Adhesive Border (with tape)
- ▶ (1) 1817x (Clamp Closure) Ultra-Clear Pouch Film
- ▶ Single-Use Package
- ▶ Clamp Included

HCPCS: A5063 + A4409

Ultra-Clear Stock No	Flange Size	For Stomas Up To
19802	1-3/4" (44 mm)	1-1/4" (32 mm)
19803	2-1/4" (57 mm)	1-3/4" (44 mm)
19804	2-3/4" (70 mm)	2-1/4" (57 mm)

## Non Sterile Kits — Urostomy



Ceramide

## 1990x

Box of 5

Urostomy Kit includes:

- ▶ (1) 1120x CeraPlus™ (extended wear) Skin Barrier
- ▶ Flat Skin Barrier with Adhesive Border (with tape)
- ▶ (1) 1892x (Ultra-Clear Pouch Film)
- ▶ (1) 7331 Urostomy Drain Tube Adapter
- ▶ Single-Use Package

HCPCS: A4432 + A4409

Ultra-Clear Stock No	Flange Size	For Stomas Up To
19902	1-3/4" (44 mm)	1-1/4" (32 mm)
19903	2-1/4" (57 mm)	1-3/4" (44 mm)
19904	2-3/4" (70 mm)	2-1/4" (57 mm)

## Hospital Kits (continued)

## One-Piece Flat Non-Sterile Kits

Box of 5 Kits

## Ultra-Clear

Kit Stock No	For Stomas Up To	Barrier Formulation	Closure Type	Pouching System Included	HCPCS	CeraPlus™ Product Alternative
89001	2-1/2" (64mm)	SoftFlex - No Tape	LNR	88401	A5061	—
89002	2-1/2" (64mm)	SoftFlex - No Tape	CLAMP	88501	A5061	—
89003	2-1/2" (64mm)	Flexend	CLAMP	8651	A4388	—
89004	2-1/2" (64mm)	Flexend	LNR	8531	A4388	—
89005	2-1/2" (64mm)	Flexend	URO	8450	A4428	—

## Two-Piece Flat Non-Sterile Kits

Box of 5 Kits

## Ultra-Clear

Kit Stock No	Flange Size	Barrier Formulation	Closure Type	Pouching System Included	HCPCS	CeraPlus™ Product Alternative
19002	1-3/4" (4mm)	Flexwear	LNR	1420x/1800x	A5063 + A4414	—
19003	2-1/4" (57mm)	Flexwear	LNR	1420x/1800x	A5063 + A4414	—
19004	2-3/4" (70mm)	Flexwear	LNR	1420x/1800x	A5063 + A4414	—
19006	4" (102mm)	Flexwear	LNR	1420x/1800x/7767*	A5063 + A4415	—
19102	1-3/4" (4mm)	Flexwear	CLAMP	1420x/1817x	A5063 + A4414	—
19103	2-1/4" (57mm)	Flexwear	CLAMP	1420x/1817x	A5063 + A4414	—
19104	2-3/4" (70mm)	Flexwear	CLAMP	1420x/1817x	A5063 + A4414	—
19106	4" (102mm)	Flexwear	CLAMP	1420x/1817x/7767*	A5063 + A4415	—
19202	1-3/4" (4mm)	Flexend	URO	1460x/1842x/7331	A4432 + A4409	19902
19203	2-1/4" (57mm)	Flexend	URO	1460x/1842x/7331	A4432 + A4409	19903
19204	2-3/4" (70mm)	Flexend	URO	1460x/1842x/7331	A4432 + A4409	19904

\*7767 is a loop ostomy bridge

## Hospital Kits (continued)

### Two-Piece Flat Sterile Kits

Box of 5 Kits

#### Ultra-Clear

Kit Stock No	For Stomas Up To	Barrier Formulation	Closure Type	Pouching System Included	HCPCS	CeraPlus™ Product Alternative
19052	1-3/4" (4mm)	Flexwear	LNR	1420x/1800x	A5063 + A4414	—
19053	2-1/4" (57mm)	Flexwear	LNR	1420x/1800x	A5063 + A4414	—
19054	2-3/4" (70mm)	Flexwear	LNR	1420x/1800x	A5063 + A4414	—
19056	4" (102mm)	Flexwear	LNR	1420x/1800x/7767*	A5063 + A4415	—
19152	1-3/4" (4mm)	Flexwear	CLAMP	1420x/1817x	A5063 + A4414	—
19153	2-1/4" (57mm)	Flexwear	CLAMP	1420x/1817x	A5063 + A4414	—
19154	2-3/4" (70mm)	Flexwear	CLAMP	1420x/1817x	A5063 + A4414	—
19156	4" (102mm)	Flexwear	CLAMP	1420x/1817x/7767*	A5063 + A4415	—
19353	2-1/4" (57mm)	Flexwear - No Tape	LNR	1520x/1800x	A5063 + A4414	—
19354	2-3/4" (70mm)	Flexwear - No Tape	LNR	1520x/1800x	A5063 + A4414	—
19252	1-3/4" (4mm)	Flexextend	URO	1460x/1842x/7331	A4432 + A4409	—
19253	2-1/4" (57mm)	Flexextend	URO	1460x/1842x/7331	A4432 + A4409	—
19254	2-3/4" (70mm)	Flexextend	URO	1460x/1842x/7331	A4432 + A4409	—
19503	2-1/4" (57mm)	Flexextend - No Tape	URO	15603/18423/7331	A4432 + A4409	—

\*7767 is a loop ostomy bridge





Pouchkins  
Pediatric  
Pouching System

# Pouchkins

**Pediatric Pouching System**

Pouchkins  
Pediatric  
Pouching System

# Pouchkins™ Pediatric Pouching System

## Premie Pouch



**3777** (CUT-TO-FIT)

**Box of 15**

- ▶ Use as a Premie Closed or Drainable Pouch
- ▶ *SoftFlex™* (standard wear) Skin Barrier — Flat
- ▶ **Without Tape Border**
- ▶ Without Starter Hole
- ▶ Without Integrated Filter
- ▶ Soft Wire Closure (if drain opening is desired)

*Fun Character Stickers Included!*

**HCPCS: A5052 + A4362**

Ultra-Clear Stock No	For Stomas Up To	Pouch Length
3777	5/8" (15 mm)	4-1/2" (11.5 mm)

## Newborn Pouch



**3778** (CUT-TO-FIT)

**Box of 15**

- ▶ Newborn Pouch; Use With or Without the Drain Valve
- ▶ *SoftFlex™* (standard wear) Skin Barrier — Flat
- ▶ **Without Tape Border**
- ▶ Without Starter Hole
- ▶ Without Integrated Filter
- ▶ Soft Wire Closure (if larger drain opening is desired)

*Fun Character Stickers Included!*

**HCPCS: A5062 + A4362**

Ultra-Clear Stock No	For Stomas Up To	Pouch Length
3778	1-3/8" x 7/8" (35 mm x 22 mm)	6-1/2" (16 cm)

## One-Piece Pouches



**3796** (CUT-TO-FIT)

**Box of 10**

- ▶ Drainable Pouch
- ▶ *SoftFlex™* (standard wear) Skin Barrier — Flat
- ▶ **Without Tape Border**
- ▶ Without Starter Hole
- ▶ Without Integrated Filter
- ▶ *Lock 'n Roll™* Microseal Closure

*Fun Character Stickers Included!*

**HCPCS: A5061**

Transparent Stock No	For Stomas Up To	Pouch Length
3796	1-1/2" (38 mm)	6-1/2" (16 cm)

## One-Piece Pouches



### 3795 (CUT-TO-FIT)

- ▶ Drainable Pouch
- ▶ *SoftFlex*™ (standard wear) Skin Barrier — Flat
- ▶ Without Starter Hole
- ▶ **Without Tape Border**
- ▶ With Integrated AF300™ Filter
- ▶ *Lock 'n Roll*™ Microseal Closure

Box of 10

*Fun Character Stickers Included!*

### HCPCS: A4424

Transparent Stock No	For Stomas Up To	Pouch Length
3795	2" (51 mm)	7" (18 cm)



### 3797 (CUT-TO-FIT)

- ▶ Urostomy Pouch
- ▶ Anti-Reflux Valve
- ▶ *SoftFlex*™ (standard wear) Skin Barrier — Flat
- ▶ Without Starter Hole
- ▶ **Without Tape Border**
- ▶ With *ComfortWear*™ Panels — Body Side Only

Box of 10  
Includes 1 adapter

*Fun Character Stickers Included!*

### HCPCS: A5071

Transparent Stock No	For Stomas Up To	Pouch Length
3797	1-1/2" (38 mm)	8-3/4" (22 cm)

## Two-Piece Skin Barrier



### 3761 (CUT-TO-FIT)

- ▶ *SoftFlex*™ (standard wear) Skin Barrier — Flat
- ▶ Without Starter Hole
- ▶ **Without Tape Border**

Box of 5

### HCPCS: A4414

Stock No	Flange Size	For Stomas Up To
3761	1-3/4" (44 mm)	1-1/4" (32 mm)

## Two-Piece Pouches



### 3799

- ▶ Drainable Pouch
- ▶ Without Integrated Filter
- ▶ *Lock 'n Roll*™ Microseal Closure
- ▶ With Belt Tabs

Box of 10

### HCPCS: A5063

Transparent Stock No	Flange Size	For Stomas Up To
3799	1-3/4" (44 mm)	6-1/2" (16 cm)



### 18292/18282

- ▶ Drainable Pouch
- ▶ With Integrated AF300™ Filter
- ▶ *Lock 'n Roll*™ Microseal Closure
- ▶ With Belt Tabs

Box of 20

### HCPCS: A4425

Transparent Stock No	Beige Stock No	Flange Size	Pouch Length
18292	18282	1-3/4" (44 mm)	7" (18 cm)



### 18382/18392

- ▶ Closed Mini-Pouch
- ▶ With Integrated AF300™ Filter
- ▶ With Belt Tabs
- ▶ With *ComfortWear*™ Panels
  - Transparent — Body Side Only
  - Beige — Both Sides

Box of 60

### HCPCS: A4419

Transparent Stock No	Beige Stock No	Flange Size	Pouch Length
18382	18392	1-3/4" (44 mm)	7" (18 cm)

# Pouchkins™ Pediatric Pouching System

## Two-Piece Pouches (continued)




**1890x/1891x**

- ▶ Urostomy Pouch
- ▶ Soft, Streamline Tap
- ▶ Drain Valve Indicator
- ▶ Belt Tabs

**Box of 10  
Includes 2 adapters**

**HCPCS: A4432**

Transparent Stock No	Beige Stock No	Flange Size	Pouch Length
18902	18912	 1-3/4" (44 mm)	9" (23 cm)

## Pediatric Ostomy Belt



**3774**

- ▶ Ostomy Belt  
- White Color

**Single Unit**

**HCPCS: A4367**

Stock No	Length
3774	10" - 17" (25 cm-43 cm) Adjustable



Accessories

Accessories

## Why use an ostomy accessory?

Ostomy accessories are used together with a pouching system. You may use one or more of these products. Here are a few reasons to use an ostomy accessory:

- It may help enhance the performance of a pouching system. For example, an accessory could help increase wear time, protect the skin, or make the removal of the adhesive easier
- It may help solve a problem such as leakage or skin irritation
- It may help improve the quality of life by eliminating odor
- Designed to meet the unique needs of your patient's body – while providing skin health, added security and confidence they deserve.

# Count on the Right Fit With CeraRing™ Barrier Rings



Hollister understands that no two people or stomas are alike — and that barrier rings can play a critical role in helping ensure the best outcome for your patients.

Whether you choose flat or convex, slim or regular, round or oval—CeraRing™ barrier rings:

- Are designed to be stretched or shaped or cut into pieces to use as a filler in uneven areas
- Help prevent stoma output from getting under the pouch seal, which may reduce leakage and skin irritation
- Can be cut or stacked together to improve the fit of the pouching system
- Are alcohol free, so there is no sting from alcohol when applied to broken or irritated skin

**Talk to your Hollister ostomy representative about the full line of CeraRing™ barrier rings with ceramide.**



## CeraRing™ Barrier Rings



Ceramide

**88xx** **Box of 10**

- ▶ Barrier Rings (Flat)
- ▶ Flexend M™ (extended wear) skin barrier with Ceramide

**HCPCS: A4385**

Stock No	Outer Diameter	Inner Diameter	Thickness
8805	2" (48 mm)	7/10" (18 mm)	4.5 mm
8815	2" (48 mm)	7/10" (18 mm)	2.3 mm <b>SLIM</b>



CONVEX  
Ceramide

**895xx** **Box of 10**

- ▶ CeraRing™ Round Convex Barrier Rings
- ▶ Flexend M™ (extended wear) skin barrier with Ceramide

**HCPCS: A4411**

Stock No	Inner Diameter
89520	13/16" (20 mm) — can be stretched to 1" (25 mm)
89530	1-3/16" (30 mm) — can be stretched to 1-3/8" (35 mm)
89540	1-9/16" (40 mm) — can be stretched to 1-3/4" (44 mm)



CONVEX  
Ceramide

**8960x** **Box of 10**

- ▶ CeraRing™ Oval Convex Barrier Rings
- ▶ Flexend M™ (extended wear) skin barrier with Ceramide

**HCPCS: A4411**

Stock No	Inner Diameter
89601	7/8" x 1-1/2" (22 mm x 38 mm) — can be stretched to 1-1/8" x 1-3/4" (27 mm x 44 mm)
89602	1-3/16" x 1-7/8" (30 mm x 48 mm) — can be stretched to 1-3/8" x 2-1/8" (35 mm x 53 mm)
89603	1-1/2" x 2-3/16" (38 mm x 56 mm) — can be stretched to 1-3/4" x 2-3/8" (44 mm x 61 mm)

## CeraPlus™ Barrier Extenders



Ceramide

**79402** **Box of 30**

- ▶ Ceramide Infused Barrier Extenders
- ▶ Hydrocolloid Formulation

**HCPCS: A4362\***

Stock No
79402

\* 3 extenders = 1 unit, Box of 30 = 10 billable units

## CeraPlus™ Skin Barrier




Ceramide

**8900** **Box of 5**

- ▶ Ceramide infused barrier formulation
- ▶ No starter hole
- ▶ No adhesive (tape)

**HCPCS: A4385**

Stock No	Size
8900	4" x 4" (10 cm x 10 cm)

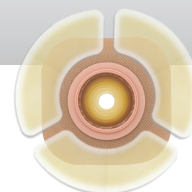

## A Sense of Security Means Everything

Three strips create a custom frame around your skin barrier.

### CeraPlus™ Barrier Extenders

At Hollister™ Ostomy Care, we know a sense of security and leakage prevention are top priorities for you.

Our CeraPlus™ barrier extenders offer an additional layer of protection for those concerned about lifting around the edges of their skin barrier. With a thin and flexible hydrocolloid formulation, CeraPlus™ barrier extenders conform easily to body contours. Three convenient strips make it easy to create a custom frame around the skin barrier. The end result is enhanced sense of security, when it matters most.

**Adapt™ Barrier Rings**



**7806** **Box of 10**

- ▶ *Flexend M™ (extended wear) skin barrier*

**HCPCS: A4385**

Stock No	Outer Diameter	Thickness
7806	4" (98 mm)	2.3 mm

**Adapt™ Skin Barrier Strips**



**79400** **Box of 10**

- ▶ Skin Barrier Strips
- ▶ *Flexend M™ (extended wear) skin barrier*

**HCPCS: A4385**

Stock No	Size	Length
79400	60 g	5-1/8" (13 cm)

**Adapt™ Skin Barrier Paste**



**79300/79301**

- ▶ Non-Adhesive
- ▶ Contains Alcohol

**HCPCS: A4406**

Stock No	Description	Quantity
79300	2.1 oz (60 g) tube	Box of 1
79301	0.5 oz (14 g) tube	Box of 20

**Adapt™ Stoma Powder**



**7906** **Box of 1**

**HCPCS: A4371**

Stock No	Description
7906	1 oz (28.3 g) puff bottle

**Adapt™ Odor Eliminator and Lubricant**



**78500/78501**

**HCPCS Code: A4394**

Stock No	Description	Quantity
78500	8 oz (236 mL) bottle	Box of 1
78501	0.27 oz (8 mL) packet	50 Packets/Box

**Adapt™ Ostomy Belts**



**7300/7299** **Box of 10**

**HCPCS: A4367**

Beige Stock No	Size
7300	23"-43" (58 cm-109 cm) Medium - Adjustable
7299	34"-65" (86 cm-165 cm) Large - Adjustable

**Adapt™ Skin Prepping Product**



**7917**

- ▶ Skin Protective Wipes
- ▶ Non-Adhesive
- ▶ Alcohol Free

**HCPCS: A5120**

Stock No
7917

**1 Box**  
**50 Wipes Each**

**NO STING**

**Adapt™ Adhesives and Removers**



**7760**

- ▶ Universal Remover Wipes
- ▶ For Removing Adhesives and Barriers
- ▶ Silicone Based
- ▶ Alcohol Free

**HCPCS: A4456**

Stock No
7760

**1 Box**  
**50 Wipes Each**

**NO STING**



**7737**

- ▶ Medical Adhesive Remover Spray
- ▶ For Removing Adhesives and Barriers
- ▶ Silicone Based
- ▶ Alcohol Free

**HCPCS: A4455**

Stock No	Description
7737	1.7 oz (50 mL) 360° spray can
7731	3.4 oz (100 mL) 360° spray can

**Box of 1**

**NO STING**

**Flexlend™ (extended wear) Skin Barrier**



**8800**  
**8801**  
**HCPCS for 8800: A4385**  
**HCPCS for 8801: A5122**

Stock No	Size	Quantity
8800	4" x 4" (10 cm x 10 cm)	Box of 5
8801	8" x 8" (20 cm x 20 cm)	Box of 3

**Skin Barrier**



**7700** **Box of 5**

▶ *HolliHesive™* (standard wear) Skin Barrier

**HCPCS: A4362**

Stock No	Size
7700	4"x 4" (10 cm x 10 cm)

**Odor Eliminators**



**773x**

- ▶ *m9™* Odor Eliminator Room Spray
- ▶ Neutralizes Odor Rather Than Masking It

**HCPCS: No Code**

Stock No (Unscented)	Description	Quantity
7732	2 oz (60 mL) pump spray	Box of 12
7733	8 oz (240 mL) pump spray	Box of 6

Stock No (Green Apple Scented)	Description	Quantity
7734	2 oz (60 mL) pump spray	Box of 12
7735	8 oz (240 mL) pump spray	Box of 6

**771x**

- ▶ *m9™* Odor Eliminator Drops

**HCPCS: A4394**



Stock No	Description	Quantity
7715	1 oz (30 mL) bottle	Box of 12
7717	8 oz (240 mL) bottle	Box of 6

**Cleaner/Decrystallizer**



**7736**

**Box of 12**

- ▶ *m9™* Cleaner/Decrystallizer
- ▶ Cleans Urinary Drainage Systems
- ▶ pH Balanced and Non-Acidic
- ▶ Free Wash Bottle Included

**HCPCS: A5131**

Stock No	Description
7736	16 oz (480 mL) bottle

**Drainable Pouch Clamp**



**8770**

**Box of 20**

**HCPCS: A4363**

Stock No	Color
8770	Beige

**Urostomy Drain Tube Adapter**



**7331**

**Box of 10**

**HCPCS: A4421**

Stock No
7331

**Hollister Universal Catheter Access Port**



**9779**

**Box of 8**

- ▶ Graduated nipple and punch allows introduction of catheter or drain through pouch
- ▶ Sterile

Stock No
9779

## Other Accessories

### Loop Ostomy Bridge



**7767**

**Box of 20**

- ▶ Sterile

**HCPCS: A4421**

**Stock No**

7767

### Bedside Drainage Collection System with Anti-Reflux Valve



**9839**

**Box of 20**

- ▶ Sample Port
- ▶ Sterile Fluid Pathway

**HCPCS: A4357**

Stock No	Description
9839	68 oz (2000 mL) Capacity

### High Output Collection Bag



**5551**

**Box of 10**

- ▶ 2L Collection Bag
- ▶ To Be Used With 1802x/830xx/26064 Pouches

**HCPCS: A4413**

**Stock No**

5551

### Vented Urinary Leg Bag Combination Pack



**9645/9655**

**Box of 10**

- ▶ Includes:
  - (1) Sterile, Vented Urinary Leg Bag — Medium or Large
  - (1) Sterile, Extension Tubing with Vacuum Relief Valve, (18" Long, Oval, 11/32" Effective Inside Diameter), with Connector
  - (1) Pair, Fabric Leg Bag Straps — Medium (23" Long, 1-1/8" Wide)
- ▶ Not Made with Natural Rubber Latex

**HCPCS: A435801**

Stock No	Description	Width
9645	Medium – 18 oz (540 mL) 10"	4-1/4"
9655	Large – 30 oz (900 mL) 11-1/8"	5"

### Urinary Leg Bag Combination Pack



**9348/9349**

**Box of 12**

- ▶ Includes:
  - (1) Sterile, Urinary leg Bag — Medium or Large
  - (1) Non-Sterile, Extension Tubing (18" Long, Oval, 11/32" Effective Inside Diameter) with Connector
  - (1) Pair, Fabric Leg Bag Straps — Medium (23" Long, 1-1/8" Wide)

**HCPCS: A4358**

Stock No	Description	Length	Width
9348	Medium — 18 oz (540 mL)	10"	4-1/4"
9349	Large — 30 oz (900 mL)	11-1/8"	5"

### Urinary Leg Bag



**9814/9805**

**Box of 10**

- ▶ Includes:
  - (1) Sterile, Urinary Leg Bag — Medium or Large
- ▶ Not Made with Natural Rubber Latex

**HCPCS: A9270**

**HCPCS: A4358 (if straps are added)**

Stock No	Description	Length	Width
9814	Medium — 18 oz (540 mL)	10"	4-1/4"
9805	Large — 30 oz (900 mL)	11-1/8"	5"

### Fabric Leg Bag Straps



**9342/9343**

**10 Pair**

- ▶ Not Made with Natural Rubber Latex

**HCPCS: A5114**

Stock No	Description	Length	Width
9842	Small	15"	1-1/8"
9343	Medium	23"	1-1/8"

### Extension Tubing



**9345/9346**

**Box of 10**


- ▶ Kink-resistant with Connector
- ▶ Oval, 11/32" Effective Inside Diameter
- ▶ Thermoplastic Elastomer
- ▶ Not Made with Natural Rubber Latex

**HCPCS: A4331**

Stock No	Description	Length
9345	Non Sterile	18"
9346	Sterile	18"

## General Information

- ▶ Similar to an enema with the water being placed into the stoma.
- ▶ Routine irrigation is recommended to be done on doctor's orders (to regulate emptying the bowel at a scheduled time).
- ▶ Appropriate for some people with sigmoid or descending colostomies.
- ▶ Not appropriate for people with urostomy or ileostomy.
- ▶ Allow approximately one hour for irrigation procedure.

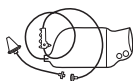

### 7721

- ▶ Cone Irrigator Kit Includes:
  - ▶ (1) Irrigator
  - ▶ (1) Stoma Cone with Connector
  - ▶ (1) Cleaning Brush

Box of 1

**Irrigation Bag Holds  
Approximately  
1500 CC**

HCPCS: A4398 + A4399

Stock No
7721



### 7723

- ▶ Stoma Cone with Connector
- ▶ Replacement Unit

Box of 10

HCPCS: A4399

Stock No
7723



### 1815x

- ▶ Drainable Irrigator Sleeve to be used with Two-Piece System
- ▶ Fits any barrier that has a 1-3/4", 2-1/4" or 2-3/4" flange
- ▶ With Belt Tabs
- ▶ Includes Clamp

Box of 5

HCPCS: A4397

Transparent Stock No	Flange Size	Sleeve Length
18152	 1-3/4" (44 mm)	30" (76 cm)
18153	 2-1/4" (57 mm)	30" (76 cm)
18154	 2-3/4" (70 mm)	30" (76 cm)



### 7728/7724

Box of 20

- ▶ Irrigator Sleeve with Belt Tabs
- ▶ Includes Clamp
- ▶ A belt is required to secure the sleeve to the body
  - \* See page 58 for belts

HCPCS: A4397

Stock No	Opening Size	Sleeve Length
7728	2" (51 mm)	22" (56 mm)
7724	3" (76 mm)	22" (56 mm)



### 1796 (PRE-SIZED)

Box of 30

- ▶ Stoma Cap
- ▶ *SoftFlex* (standard wear) Skin Barrier — Flat
- ▶ Without Tape Border
- ▶ With Integrated Filter

HCPCS: A5055

Beige Stock No	Barrier Opening	Size
1796	1-15/16" (50 mm)	4" (102 mm)



### 318x (PRE-SIZED)

Box of 30

- ▶ Stoma Cap
- ▶ Adhesive Border (with tape)
- ▶ With Integrated Filter
- ▶ Odor-Barrier Film

HCPCS: A5055

Beige Stock No	Opening	Size
3184	2" (51 mm)	4-1/4" (11 cm)
3186	3" (76 mm)	4-1/4" (11 cm)







Education

Education









### Ostomy Clinical Education

Knowledge is power. When you have access to the latest research, educational materials, product evidence, and training you can do your best work. And that gives your patients the best chance for a successful ostomy journey.

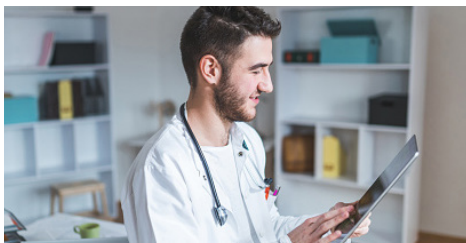


### Continued Learning

Ostomy care is constantly changing and evolving, requiring you to be a student throughout your career. We are committed to helping you stay current. Here, you can access webinars, online courses, case studies, and other tools that can help you acquire knowledge and learn new skills – at your own pace.

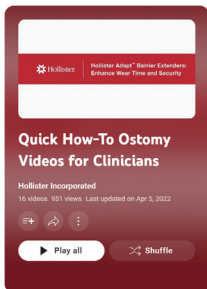

### Quick How-To Ostomy Videos for Clinicians

Provides brief instructions for use of specific products and/or instructional information related to specific ostomy products or procedures.

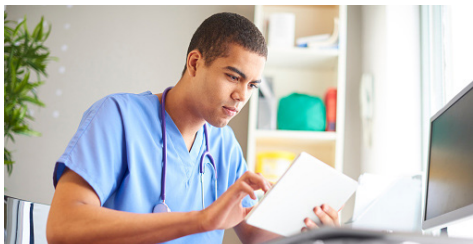

### Quick Educational Tools

Access helpful videos, guides, and care tips to learn more about what to expect with ostomy surgery, how to use ostomy products, maintain healthy skin around the stoma, and resume lifestyle activities with an ostomy.




### Ostomy Care Tips

These ostomy care tips are designed to assist with ostomy care education. Topics include general ostomy care, product and accessories instructions, maintaining healthy peristomal skin, insurance information, and more.



### International Resources

We have translated our educational resources into a variety of languages. These materials can help your patients around the world understand their type of ostomy and how they can live their lives to the fullest.

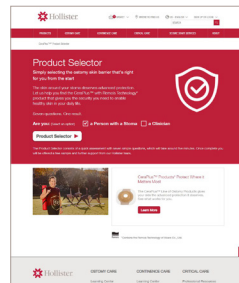


### Product Selector Guide

By answering seven simple questions, let us help you find the CeraPlus™ Ostomy Product\* you need. Once completed, you will be offered a free sample and further support from our Hollister Team.



### Peristomal Skin Assessment Guide

The Peristomal Skin Assessment Guide helps with assessment and management of common peristomal skin complications and is available on the WOCN Society website.

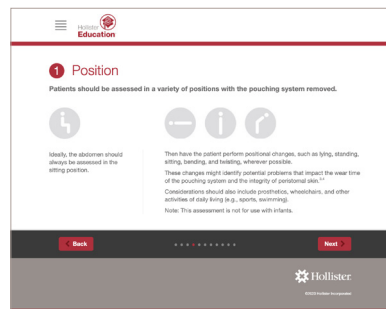

### Fit Indicator Tool (FIT)

Fit Indicator Tool (FIT) has been designed to help provide a simple patient assessment methodology to determine the need for a flat or convex barrier.

The following items can be ordered through your local Hollister Sales Representative:

**Stoma Measuring Guides**

Stoma Measuring Guide 903515  
Oval Stoma Measuring Guide 923298  
(Printable versions available online at

[www.hollister.com/en/ostomycare/educationaltools](http://www.hollister.com/en/ostomycare/educationaltools)

**Pre-Surgery Care Kits**

Includes guidebook to facilitate the patient's pre-surgery visit with a clinician and materials for self-led practice applying and removing an ostomy pouching system with a model stoma.

Colostomy Pre-Surgery Care Kit K2025

Ileostomy Pre-Surgery Care Kit K2026

Urostomy Pre-Surgery Care Kit K2027

**Patient Discharge Folders**

Colostomy Patient Discharge Folder 923261

Ileostomy Patient Discharge Folder 923262

Urostomy Patient Discharge Folder 923263

Scan here for additional Hollister product offerings



Ostomy Care

---



Continence Care

---



Critical Care

---

<u>Stock No</u>	<u>Page(s)</u>	<u>Stock No</u>	<u>Page(s)</u>	<u>Stock No</u>	<u>Page(s)</u>	<u>Stock No</u>	<u>Page(s)</u>	<u>Stock No</u>	<u>Page(s)</u>	<u>Stock No</u>	<u>Page(s)</u>
1120x	16	1801x	45	1925x	50	7736	61	828x	39	8800	60
11201x	16	1802x	45	1935x	50	7737	59	829x	39	8801	60
1140x	15	1810x	19	1950x	50	7760	59	830xx	46	883xx	42
1150x	15	1811x	18	1960x	21, 48	7767	62	8331	40	884xx	42
1170x	15	1812x	19	1980x	21, 48	7806	58	8371	40	88501	42
1190x	15	1813x	18	1990x	21, 48	78500	58	8381	41	887xx	42
1340x	36	1814x	18	26064	46	78501	58	8413xx	27	88xx	56
1350x	36	1815x	64	318x	65	7906	58	841x	27	88800	42
139xx	34	1816x	18	3761	53	7917	59	8429x	28	8900	56
1420x	36	1817x	19	3774	54	79300	58	843911	38	8900x	49
1430x	36	1818x	18	3777	51	79301	58	8440xx	28	8901	25
1460x	36	1819x	18	3778	51	79400	58	8449x	40	891xx	24
1470x	35	1820x	19	3795	52	79402	56	84590	40	89301	26
1480x	34	1828x	18, 53	3796	51	80070	46	8469x	40	8931	26
1510x	16	1829x	18, 53	3797	52	80110	46	847xx	38	89311	25
15111x	16	1835x	20	3799	53	8031	41	8479xx	38	89321	25
1520x	36	1836x	19	5551	62	8081	41	8489x	39	895xx	23, 56
1530x	15	1837x	19	7299	58	8181	41	8499x	39	89518	23
15310x	16	1838x	19	7300	58	8195111	23	8531	40	8960x	56
1540x	36	1839x	19	7331	61	8195411	23	854x	38	896x	24
1550x	36	1875x	20	7700	60	821xx	42	8551	40	896xx	23
1560x	35	1890x	20, 54	771x	61	823xx	42	855x	40	898xx	24
1580x	34	1891x	20, 54	7721	64	824x	38	857xx	38	89xx	26
1590x	34	1900x	49	7723	64	824xx	42	863x	40	899x	25
1610x	35	1905x	50	7724	65	82500	42	864x	40	977x	47
1640x	36	1910x	49	7728	65	8261x	39	866x	38	9779	61
1796	65	1915x	50	773x	61	8262x	38	867xx	38	9839	62
1800x	18	1920x	49	7731	59	827xx	38	8770	61		

## HCPCS Codes Disclaimer

The reimbursement information provided herein is intended to provide general information concerning coding of Hollister products only. Hollister does not guarantee coverage or payment for any products. The ultimate responsibility for proper coding, satisfying reimbursement requirements, and obtaining reimbursement remains with the provider. Coding and coverage policies and guidelines are complex, can vary from one carrier or region to another, and are updated frequently. Providers should check with their local carriers or intermediaries often and should consult with counsel, a reimbursement specialist, and/or the DMECS website for any coding, coverage, reimbursement, or billing questions.

*HCPCS codes shown in this Product Catalog are current as of January 2024.*

## General Medicare Product Maximum Allowables of Selected Ostomy Products

Two-Piece Pouching Systems	Maximum Quantity per Month
Flat, Standard Wear Skin Barriers with Floating Flanges	20
Drainable Pouches	20
Closed Pouches	60
Urostomy Pouches	20
Irrigator Sleeves	4

One-Piece Pouching Systems	Maximum Quantity per Month
Drainable Pouches with Flat, Standard Wear Skin Barriers	20
Closed Pouches with Flat, Standard Wear Skin Barriers	60
Stoma Caps	31

Pediatric Pouching Systems	Maximum Quantity per Month
One-Piece Drainable and Urostomy Pouches	20
Two-Piece Skin Barriers	20
Two-Piece Drainable and Urostomy Pouches	20
Two-Piece Mini Closed Pouches	60
Pediatric Ostomy Belt	1

Accessories	Maximum Quantity per Month
Skin Barrier Rings	*
Skin Barrier Paste	4 oz
Stoma Powder	10 oz/6 months
Lubricating Deodorant	*
Ostomy Belt	1
Skin Protective Wipes	150/6 months

**Note:** Medicare allows a usual maximum quantity of supplies each month. The table above provides the maximum allowable amounts. The amount of supplies you need may vary from the maximum quantities. That will be determined by your physician and must be on your order and in your medical record. If you require more ostomy supplies than what is listed in the table, your physician will need to provide a medical justification for the additional supplies. The ostomy supplier will need to have that medical justification on file when billing the Medicare program for the additional supplies.

\* No usual maximum amount established.

## Reimbursement For Fistula Management by Medicare

Reimbursement for fistula management under the Medicare program is described in the Surgical Dressings Policy Article. Specifically the article states:

“Surgical dressings are covered when either of the following criteria are met:

1. They are required for the treatment of a wound caused by, or treated by, a surgical procedure; **or**
2. They are required after debridement of a wound.”

The article goes on to explain:

“Examples of situations in which dressings are noncovered under the Surgical Dressings benefit are:

- a. Drainage from a cutaneous fistula which has not been caused by or treated by a surgical procedure”

For information see The Fistula Management Care Tip on [www.hollister.com](http://www.hollister.com)

## Latex Statement

All Hollister Ostomy Care products found in this catalog and its packaging does not include natural rubber latex as a component, and our component suppliers do not add natural rubber as part of their production process.

Prior to use, be sure to read the Instructions for Use for information regarding intended use, contraindications, warnings, precautions and instructions.

## Medicare's Usual Maximum Quantity for Ostomy Supplies

The quantity of ostomy supplies needed by a patient is determined primarily by the type of ostomy, its location, its construction, and the condition of the skin surface surrounding the stoma. There will be variation according to individual patient need and their needs may vary over time. The table below lists the maximum number of items/units of service that are usually medically necessary. The actual quantity needed for a particular patient may be more or less than the amount listed depending on the factors that affect the frequency of barrier and pouch change. The medical necessity for use of a greater quantity of supplies than the amounts listed must be clearly documented in the patient's medical record, and must be available upon request. If adequate documentation is not provided when requested, the excess quantities will be denied as not medically necessary.

**Per Medicare LCD for Ostomy Supplies (L33828), effective 1-1-2024.**

HCPCS	Long Description	Maximum Quantity
A4331	Extension drainage tubing, any type, any length, with connector/adaptor, for use with urinary leg bag or urostomy pouch, each	*
A4357	Bedside drainage bag, day or night, with or without anti-reflux device, with or without tube, each	2 ea /mo
A4361	Ostomy faceplate, each	3 ea/6 mo
A4362	Skin barrier, solid, 4 x 4 or equivalent, each	20 ea/mo
A4363	Ostomy clamp, any type, replacement only, each	*
A4364	Adhesive, liquid or equal, any type, per ounce	4 oz/mo
A4366	Ostomy vent, any type, each	*
A4367	Ostomy belt, each	1 ea/mo
A4368	Ostomy filter, any type, each	*
A4369	Ostomy skin barrier, liquid (spray, brush, etc.), per ounce	2 oz/mo
A4371	Ostomy skin barrier, powder, per ounce	10 oz/6 mo
A4372	Ostomy skin barrier, solid 4 x 4 or equivalent, standard wear with built-in convexity, each	*
A4373	Ostomy skin barrier, with flange (solid, flexible, or accordion), with built-in convexity, any size, each	*
A4375	Ostomy pouch, drainable, with faceplate attached, plastic, each	*
A4376	Ostomy pouch, drainable, with faceplate attached, rubber, each	*
A4377	Ostomy pouch drainable for use on faceplate, plastic, each	10 ea/mo
A4378	Ostomy pouch, drainable, for use on faceplate, rubber, each	*
A4379	Ostomy pouch, urinary, with faceplate attached, plastic, each	*
A4380	Ostomy pouch, urinary, with faceplate attached, rubber, each	*
A4381	Ostomy pouch, urinary, for use with faceplate, plastic, each	10 ea/mo
A4382	Ostomy pouch, urinary, for use on faceplate, heavy plastic, each	*
A4383	Ostomy pouch, urinary, for use on faceplate, rubber, each	*
A4384	Ostomy faceplate equivalent, silicone ring, each	*
A4385	Ostomy skin barrier, solid 4 x 4 or equivalent, extended wear, without built-in convexity, each	*
A4387	Ostomy pouch, closed, with barrier attached, with built-in convexity (one-piece), each	*
A4388	Ostomy pouch, drainable, with extended wear barrier attached, (one-piece), each	*
A4389	Ostomy pouch, drainable, with barrier attached, with built-in convexity (one-piece), each	*

\*No established usual maximum, refer to Medicare Utilization Guidelines.

HCPCS	Long Description	Quantity
A4390	Ostomy pouch, drainable, with extended wear barrier attached, with built-in convexity (one-piece), each	*
A4391	Ostomy pouch, urinary, with extended wear barrier attached (one-piece), each	*
A4392	Ostomy pouch, urinary, with standard wear barrier attached, with built-in convexity (one-piece), each	*
A4393	Ostomy pouch, urinary, with extended wear barrier attached, with built-in convexity (one-piece), each	*
A4394	Ostomy deodorant, with or without lubricant, for use in ostomy pouch, liquid, per fluid ounce	*
A4395	Ostomy deodorant for use in ostomy pouch, solid, per tablet	*
A4396	Ostomy belt with peristomal hernia support	*
A4397	Irrigation supply, sleeve, each	4 ea/mo
A4398	Ostomy irrigation supply, bag, each	2 ea/6 mo
A4399	Ostomy irrigation supply, cone/catheter, with or without brush	2 ea/6 mo
A4402	Lubricant, per ounce	4 oz/mo
A4404	Ostomy ring, each	10 ea/mo
A4405	Ostomy skin barrier, non-pectin based, paste, per ounce	4 oz/mo
A4406	Ostomy skin barrier, pectin-based, paste, per ounce	4 oz/mo
A4407	Ostomy skin barrier, with flange (solid, flexible, or accordion), extended wear, with built-in convexity, 4 x 4 inches or smaller, each	*
A4408	Ostomy skin barrier, with flange (solid, flexible, or accordion), extended wear, with built-in convexity, larger than 4 x 4 inches, each	*
A4409	Ostomy skin barrier, with flange (solid, flexible, or accordion), extended wear, without built-in convexity, 4 x 4 inches or smaller, each	*
A4410	Ostomy skin barrier, with flange (solid, flexible, or accordion), extended wear, without built-in convexity, larger than 4 x 4 inches, each	*
A4411	Ostomy skin barrier, solid 4 x 4 or equivalent, extended wear, with built-in convexity, each	*
A4412	Ostomy pouch, drainable, high output, for use on a barrier with flange (two-piece system), without filter, each	*
A4413	Ostomy pouch, drainable, high output, for use on a barrier with flange (two-piece system), with filter, each	*
A4414	Ostomy skin barrier, with flange (solid, flexible, or accordion), without built-in convexity, 4 x 4 inches or smaller, each	20 ea/mo
A4415	Ostomy skin barrier, with flange (solid, flexible, or accordion), without built-in convexity, larger than 4 x 4 inches, each	20 ea/mo
A4416	Ostomy pouch, closed, with barrier attached, with filter (one-piece), each	60 ea/mo
A4417	Ostomy pouch, closed, with barrier attached, with built-in convexity, with filter (one-piece), each	60 ea/mo
A4418	Ostomy pouch, closed, without barrier attached, with filter (one-piece), each	60 ea/mo
A4419	Ostomy pouch, closed, for use on barrier with non-locking flange, with filter (two-piece), each	60 ea/mo
A4420	Ostomy pouch, closed, for use on barrier with locking flange (two-piece), each	60 ea/mo
A4421	Ostomy supply; miscellaneous	*
A4422	Ostomy absorbent material (sheet/pad/crystal packet) for use in ostomy pouch to thicken liquid stomal output, each	*

\*No established usual maximum, refer to Medicare Utilization Guidelines.

# Medicare Utilization Chart

HCPCS	Long Description	Quantity
A4423	Ostomy pouch, closed, for use on barrier with locking flange, with filter (two-piece), each	60 ea/mo
A4424	Ostomy pouch, drainable, with barrier attached, with filter (one-piece), each	20 ea/mo
A4425	Ostomy pouch, drainable, for use on barrier with non-locking flange, with filter (two-piece system), each	20 ea/mo
A4426	Ostomy pouch, drainable, for use on barrier with locking flange (two-piece system), each	20 ea/mo
A4427	Ostomy pouch, drainable, for use on barrier with locking flange, with filter (two-piece system) each	20 ea/mo
A4428	Ostomy pouch, urinary, with extended wear barrier attached, with faucet-type tap with valve (one-piece), each	*
A4429	Ostomy pouch, urinary, with barrier attached, with built-in convexity, with faucet-type tap with valve (one-piece), each	20 ea/mo
A4430	Ostomy pouch, urinary, with extended wear barrier attached, with built-in convexity, with faucet-type tap with valve (one-piece), each	*
A4431	Ostomy pouch, urinary, with barrier attached, with faucet-type tap, with valve (one-piece), each	20 ea/mo
A4432	Ostomy pouch, urinary, for use on barrier with non-locking flange, with faucet-type tap with valve (two-piece), each	20 ea/mo
A4433	Ostomy pouch, urinary, for use on barrier with locking flange (two-piece), each	20 ea/mo
A4434	Ostomy pouch, urinary, for use on barrier with locking flange, with faucet-type tap with valve (two-piece), each	20 ea/mo
A4435	Ostomy pouch, drainable, high output, with extended wear barrier (one-piece system), with or without filter, each	*
A4450	Tape, non-waterproof, per 18 square inches	40 ea (18 sq in)/mo
A4452	Tape, waterproof, per 18 square inches	40 ea (18 sq in)/mo
A4455	Adhesive remover or solvent (for tape, cement, or other adhesive), per ounce	16 oz/6 mo
A4456	Adhesive remover, wipes, any type, each	*
A5051	Ostomy pouch, closed, with barrier attached (one-piece), each	60 ea/mo
A5052	Ostomy pouch, closed, without barrier attached, (one-piece), each	60 ea/mo
A5053	Ostomy pouch, closed, for use on faceplate, each	60 ea/mo
A5054	Ostomy pouch, closed, for use on barrier with flange (two-piece), each	60 ea/mo
A5055	Stoma Cap	31 ea/mo
A5056	Ostomy pouch, drainable, with extended wear barrier attached, with filter, (one-piece), each	40 ea/mo
A5057	Ostomy pouch, drainable, with extended wear barrier attached, with built-in convexity, with filter, (one-piece), each	40 ea/mo
A5061	Ostomy pouch, drainable, with barrier attached (one-piece), each	20 ea/mo
A5062	Ostomy pouch, drainable, without barrier attached (one-piece), each	20 ea/mo
A5063	Ostomy pouch, drainable, for use on barrier with flange (two-piece system), each	20 ea/mo
A5071	Ostomy pouch, urinary, with barrier attached (one-piece), each	20 ea/mo
A5072	Ostomy pouch, urinary, without barrier attached (one-piece), each	20 ea/mo
A5073	Ostomy pouch, urinary, for use on barrier with flange (two-piece), each	20 ea/mo

\*No established usual maximum, refer to Medicare Utilization Guidelines.

HCPCS	Long Description	Quantity
A5081	Continent device, plug for continent stoma	31 ea/mo
A5082	Continent device, catheter for continent stoma	1 ea/mo
A5083	Continent device, stoma absorptive cover for continent stoma	150 ea/mo
A5093	Ostomy accessory, convex insert	10 ea/mo
A5102	Bedside drainage bottle with or without tubing, rigid or expandable, each	2 ea/6 mo
A5120	Skin barrier, wipes or swabs, each	150 ea/6 mo
A5121	Skin barrier, solid, 6 x 6 or equivalent, each	20 ea/mo
A5122	Skin barrier, solid, 8 x 8 or equivalent, each	20 ea/mo
A5126	Adhesive or non-adhesive, disk or foam pad	20 ea/mo
A5131	Appliance cleaner, incontinence and ostomy appliances, per 16 ounce	1 ea (16 oz btl)/mo
A9270	Non-covered item or service	*

\*No established usual maximum, refer to Medicare Utilization Guidelines.



## DID YOU KNOW...

### Supplier Matching

Hollister Secure Start<sup>SM</sup> services can connect patients with a supplier of their choice based on their insurance plan for the products you recommend. The amount of supplies they need may vary from the maximum allowable quantities, but this is determined by a physician and should be included on the patient's order and medical record. For additional assistance, visit [www.securestartservices.com](http://www.securestartservices.com) or call us at 1.888.808.7456.





## Hollister CeraPlus™

Ostomy Products\* have received official dermatological accreditation from the **Skin Health Alliance**




## Who is the Skin Health Alliance?

The **Skin Health Alliance** is a leading international skin health organization that works closely with dermatologists, scientists and the skin care industry with a common interest in the promotion of healthy skin. The **Skin Health Alliance** awards professional accreditation and endorsement to companies, services and brands seeking independent dermatological recognition for their product research.

- **Hollister CeraPlus™** Ostomy Products\* have received official dermatological accreditation from the **Skin Health Alliance**
- In addition to the current rigorous safety testing our products go through, the SHA accreditation signifies that the **Hollister CeraPlus™** Ostomy Products have been independently verified by a specialist panel
- **Hollister CeraPlus™** Ostomy Products may help prevent skin problems and help keep skin healthy from Day 1
- By choosing **Hollister CeraPlus™** Ostomy Products, people using stoma pouches can have the confidence and pride of knowing their barrier helps protect their skin, empowering them to live their life to the fullest

CeraPlus Skin Barriers contain the Remois Technology of Alcare Co., Ltd.





The spirit of Hollister  
remains the same  
today as it did in 1921.

## Our Mission

To make life more rewarding  
and dignified for the people  
who use our products  
and services.

## Our Vision

To grow and prosper as an  
independent, employee-owned  
company and in the process  
become better human beings.

# Our Heritage

The roots of Hollister Incorporated date back to 1921, when 23-year-old entrepreneur John Dickinson Schneider founded a small printing company in Chicago, which he named JDS Printer Craftsman. From the very beginning, John Schneider's business philosophy was "only first class is good enough," and the company became known for its quality products, quality service, and quality people.

Hollister, the Hollister logo, Adapt, Adapt CeraRing, AF300, CeraPlus, ComfortWear, Flexextend, Flexextend M, FlexWear, Lock n' Roll, m9, New Image, Pouchkins, Premier, Secure Start and SoftFlex are trademarks of Hollister Incorporated.

All other trademarks and copyrights are the property of their respective owners.

© 2024 Hollister Incorporated

## Contact Us

**Monday - Friday: 8:00am - 5:00pm (CT)**

For consumer or healthcare professional questions:

**1.888.808.7456**

For dealer, distributor, institution and pharmacy questions or orders:

**1.800.323.4060**

**FAX 1.877.321.7805**

## Hollister Secure Start<sup>SM</sup> Services

### Enroll on the Go Clinician Assisted Enrollment

Access to personalized, no-cost support for ostomy patients is easier—with mobile device enrollment!  
Call us with any questions at 1.888.808.7456.



**Hollister Incorporated**  
2000 Hollister Drive  
Libertyville, Illinois 60048  
1.800.323.4060



[www.hollister.com](http://www.hollister.com)